

Online Bracket Library


Since Release 3.2.100 (233), licensed users of program versions


- OnyxCeph³™ 3D PRO
- OnyxCeph³™ LAB
- OnyxCeph³™ OMS



can synchronize and complete their local bracket table by accessing and downloading from the Onyx online bracket library.

Open Library

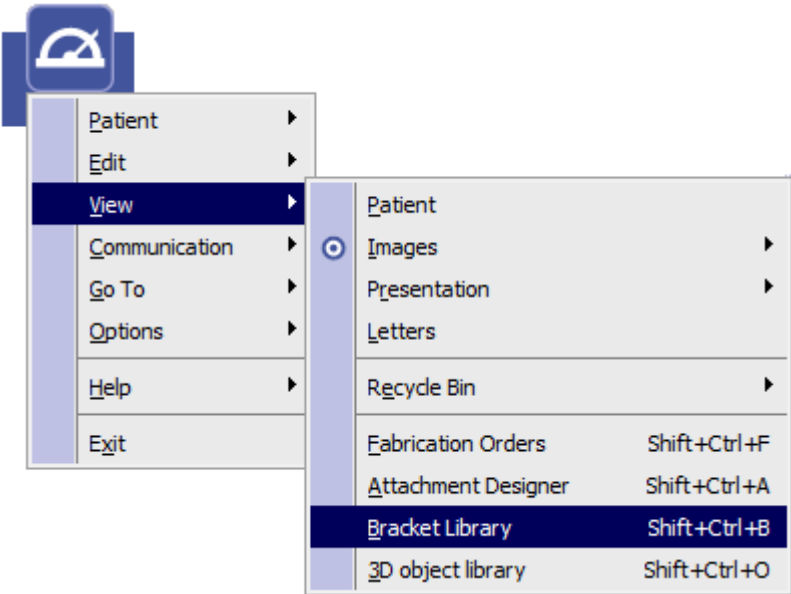The interface is available by main menu [View|BracketLibrary] or Shortcut SHIFT+CTRL+B

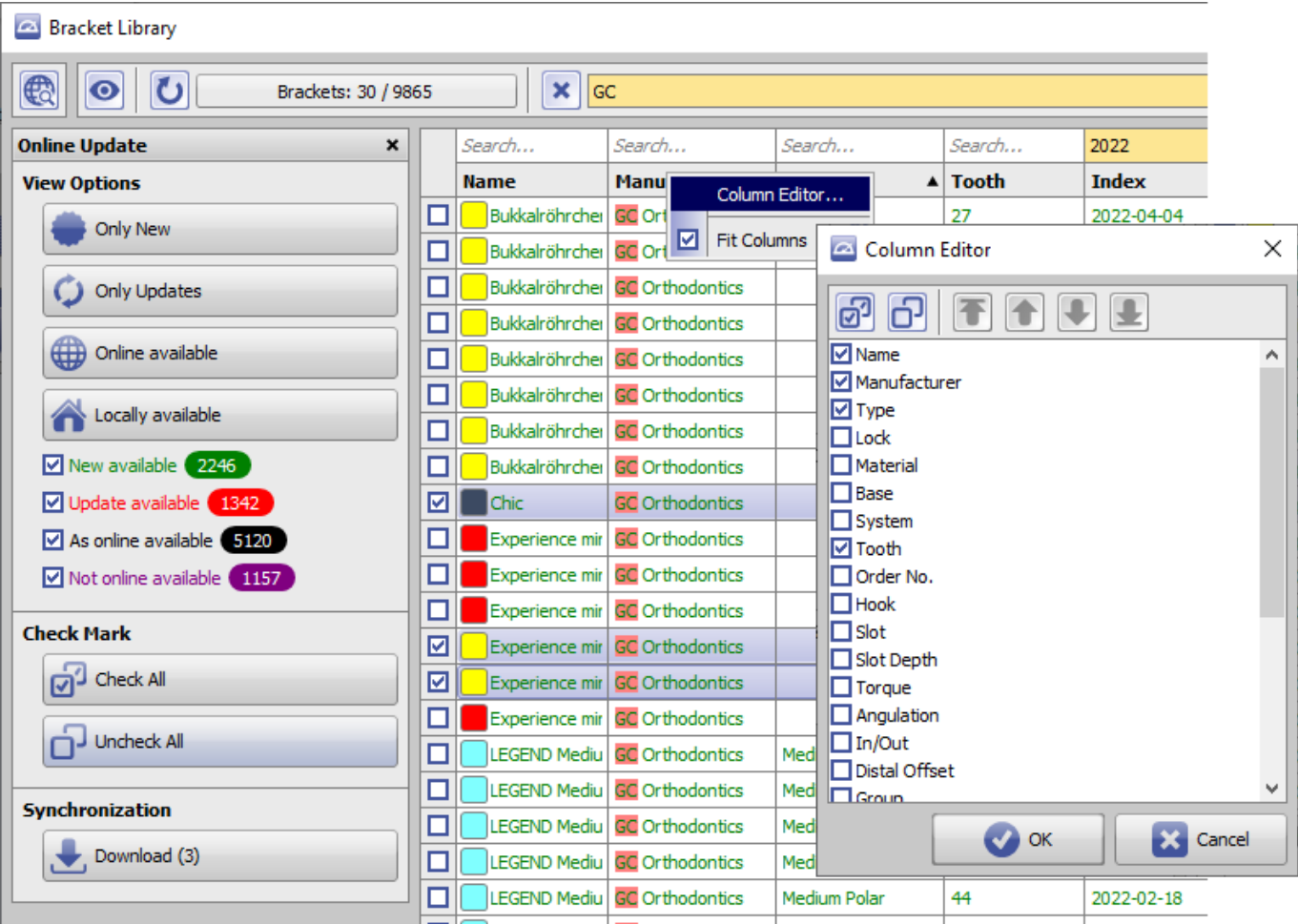



The local bracket table can be refreshed by icon button [Reload] . By icon button [Search for updates]  the content of the online bracket library can be compared with the local bracket table by the brackets ID. As result of the comparison, all bracket files available online and local are categorized and then color-coded accordingly:

Color Code	Where Is The Bracket File Located ?	Online	Local	Available For Download ?
2246	only online	✓	✗	✓ Yes
1342	online + local, different version	✓	✓	✓ Yes
5120	online + local, identical version	✓	✓	No
1157	only local	✗	✓	No

Search for Items

The bracket table display can be filtered by buttons and checkboxes in panel [View Options].

View Options

Only New

Only Updates

Online available

Locally available

☒ New available

2246

☒ Update available

1342

☒ As online available

5120

☒ Not online available

1157

Additionally, the search boxes top of each column and the full-text-search field above the table can be used to find and narrow down the list view to the searched bracket files. By column header context menu, the columns to be displayed can be configured.

×

GC

	Search...	Search...	Search...	Search...	2022
	Name	Manu		▲ Tooth	Index
<input type="checkbox"/>	Bukkalröhrchei	GC Ort		27	2022-04-04
<input type="checkbox"/>	Bukkalröhrchei	GC Ort			
<input type="checkbox"/>	Bukkalröhrchei	GC Orthodontics			
<input type="checkbox"/>	Bukkalröhrchei	GC Orthodontics			
<input type="checkbox"/>	Bukkalröhrchei	GC Orthodontics			
<input type="checkbox"/>	Bukkalröhrchei	GC Orthodontics			
<input type="checkbox"/>	Bukkalröhrchei	GC Orthodontics			

Column Editor...

☒ Fit Columns

Column Editor

☒ Name

☒ Manufacturer

☒ Type

☐ Lock

☐ Material

☐ Base

☐ System

☒ Tooth

☐ Order No.

☐ Hook

☐ Slot

☐ Slot Depth

☐ Torque

☐ Angulation

☐ In/Out

☐ Distal Offset

☐ Group

✓


OK

×

Cancel

Download Items

To finally download the online available desired brackets (marked either green or red) to the local bracket table, the checkbox left in each line needs to be activated. By buttons [Check All] or [Uncheck All] this can be done for all brackets currently displayed in the list at once:

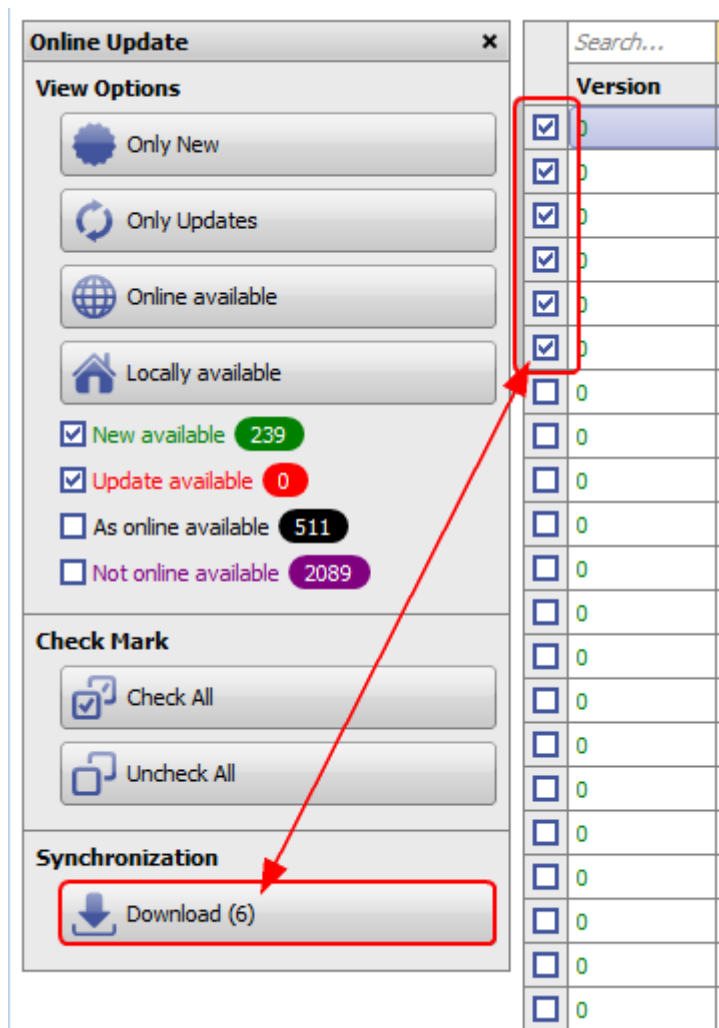

Links

- W Downloading Brackets / Attachments
- W Import Requirements

From:
<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
<http://onyxwiki.net/doku.php?id=en:3948&rev=1689140670>

Last update: **2023/07/12 07:44**

