

Module Cast Adjust



<WRAP group> <WRAP threequarters column>

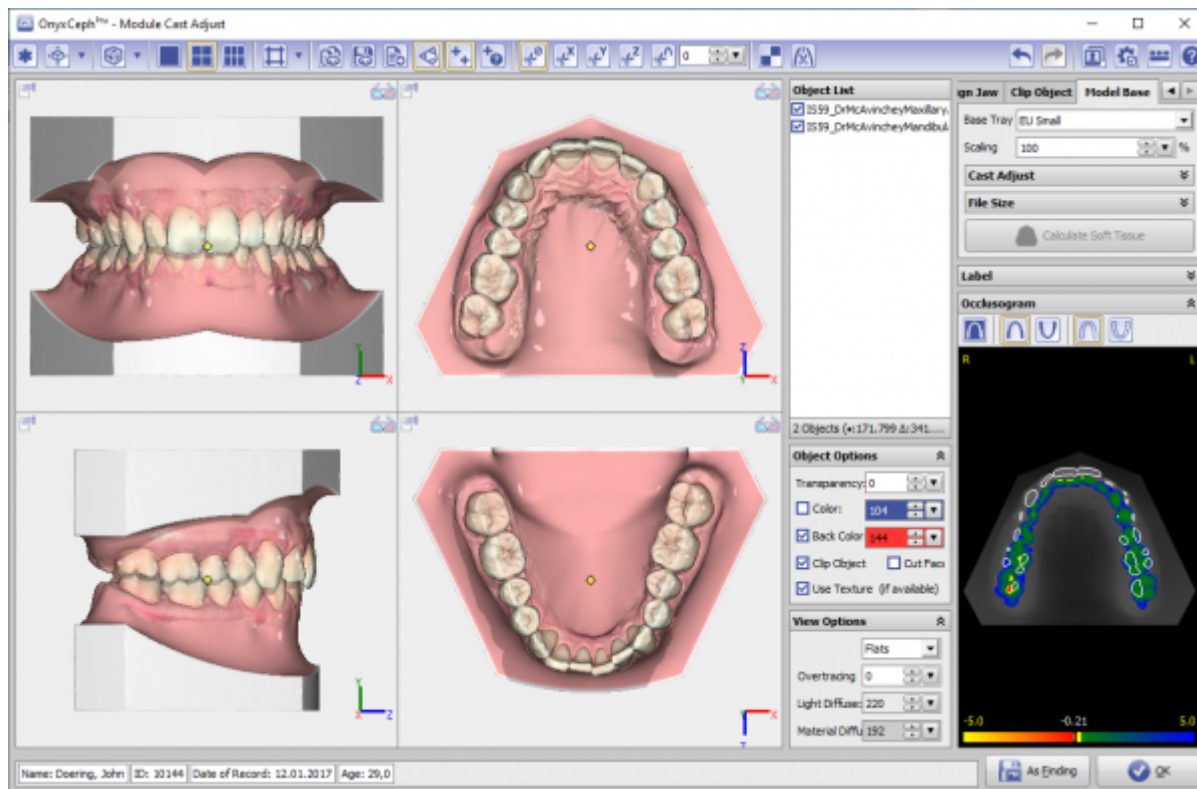
Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability		✗			

Module Cast Adjust allows to re-align the jaw relation of a scan according to different specifications and the attachment of different types of model bases (EU, US, plate, horseshoe bases) after manual definition of a trim line. The study model generated in this way must be saved as a separate finding.




First Steps

1. Open Module with scan (open shell - two arches in occlusion)
2. Optional: Adjust jaw relation on tab |Articulate|,|Jaw Rotation Axis| or |Align Jaw|
3. On tab |Clip Object| by right clicks, define trim line for upper and lower jaw
4. Clip along trim line Press
5. On tab |Model Base| select type of base tray and merge with scan
6. Save study model as new finding

Module user interface



Links

	Module Handout Cast Adjust
	https://youtu.be/TAAGkmFIKD4
	https://drive.google.com/open?id=0B5nAIDRchl5WakJBYXJ2QWNFSXM

From:
<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
http://onyxwiki.net/doku.php?id=en:base_tray&rev=1644236283

Last update: **2022/02/07 13:18**

