

# Module Cast Adjust



Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability					

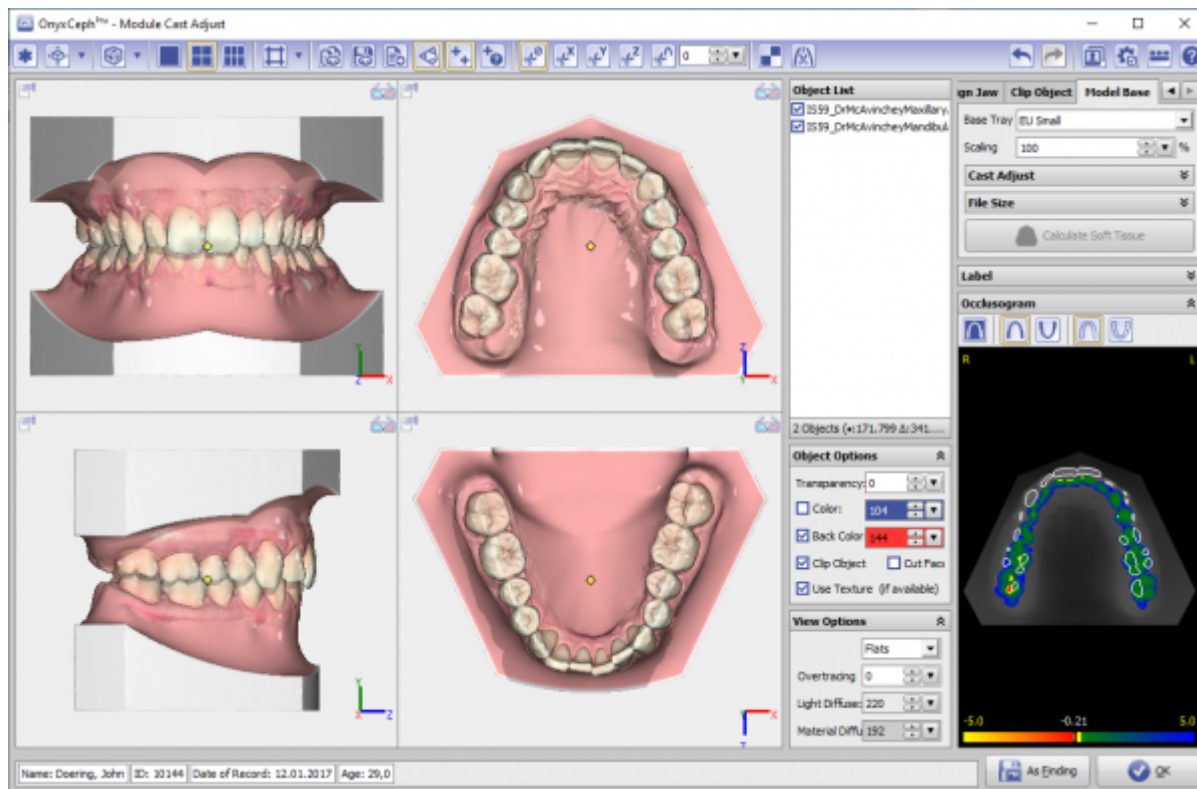
- nicht  
enthalten
- immer  
enthalten
- optional  
enthalten

Module Cast Adjust allows to re-align the jaw relation of a scan according to different specifications and the attachment of different types of model bases (EU, US, plate, horseshoe bases) after manual definition of a trim line. The study model generated in this way must be saved as a separate finding.




## First Steps

1. Open Module with scan (open shell - two arches in occlusion)
2. Optional: Adjust jaw relation on tab |Articulate|,|Jaw Rotation Axis| or |Align Jaw|
3. On tab |Clip Object| by right clicks, define trim line for upper and lower jaw
4. Clip along trim line Press
5. On tab |Model Base| select type of base tray and merge with scan
6. Save study model as new finding

## Module user interface



## Links

	<a href="#">Module Handout Cast Adjust</a>
	<a href="https://youtu.be/TAAGkmFIKD4">https://youtu.be/TAAGkmFIKD4</a>
	<a href="https://drive.google.com/open?id=0B5nAIDRchl5WakJBYXJ2QWNFSXM">https://drive.google.com/open?id=0B5nAIDRchl5WakJBYXJ2QWNFSXM</a>

From:  
<http://onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:  
[http://onyxwiki.net/doku.php?id=en:base\\_tray&rev=1644236612](http://onyxwiki.net/doku.php?id=en:base_tray&rev=1644236612)

Last update: **2022/02/07 13:23**

