

Module Bonding Jigs 3D




Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	×	×			

In module Bonding Jigs 3D module, occlusal transfer jigs can be constructed to transference the virtually planned bracket positions using metallic standard bracket adapters.

First Steps

1. Open Module with bracket-attached Model
2. Select Crowns
3. Assign adapter and Jig parameters to crowns
4. Construct Jigs individually or connected
5. Arrange Jigs for manufacturing / export
6. Save planning as project and/or as finding

Module user interface

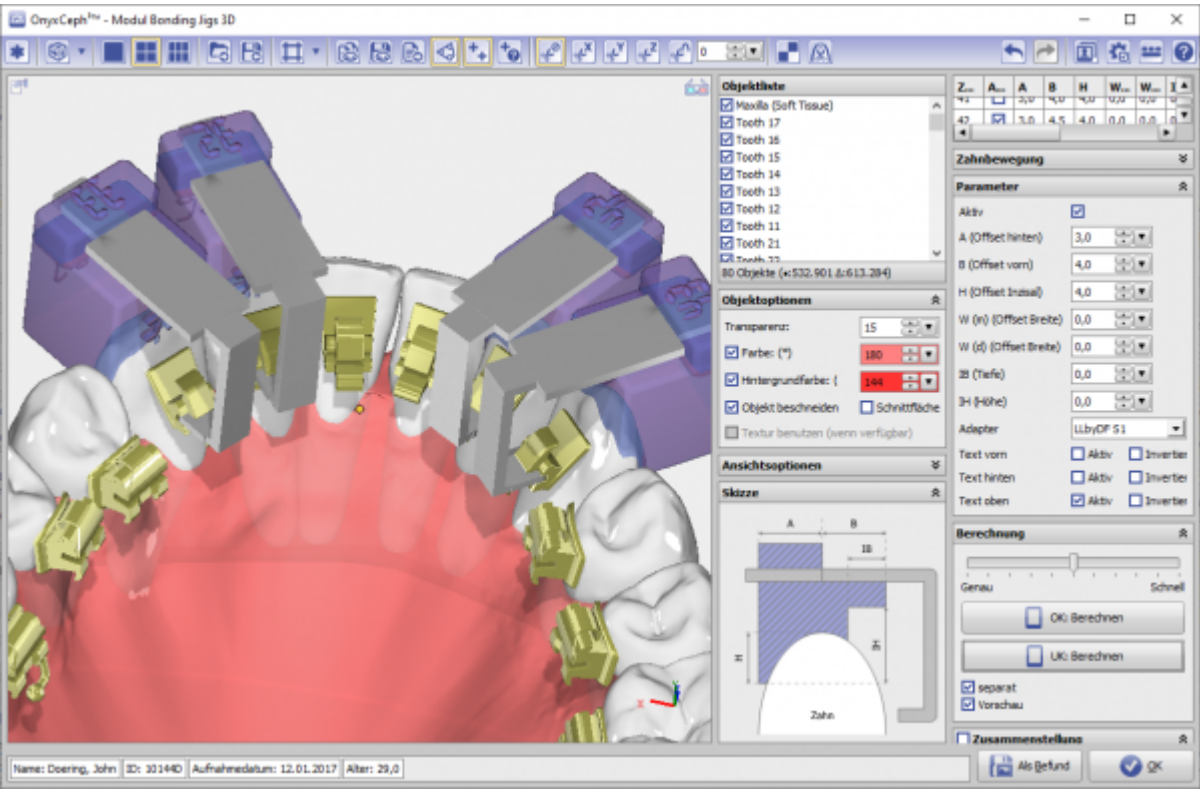

Links



From:

<http://onyxwiki.net/> - **[OnyxCeph³™ Wiki]**

Permanent link:

http://onyxwiki.net/doku.php?id=en:bonding_jigs_3d&rev=1586016179

Last update: **2020/04/04 18:02**

