

Module Digitize 2D



Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	✗	✗			

not included

included by default

optionally included

Modul Digitize 2D allows the cephalometric evaluation of 2D image data records. Depending on the image type and the type of tracing assigned to it, various pre-installed methods (analyzes) can be selected and evaluated individually or in groups.

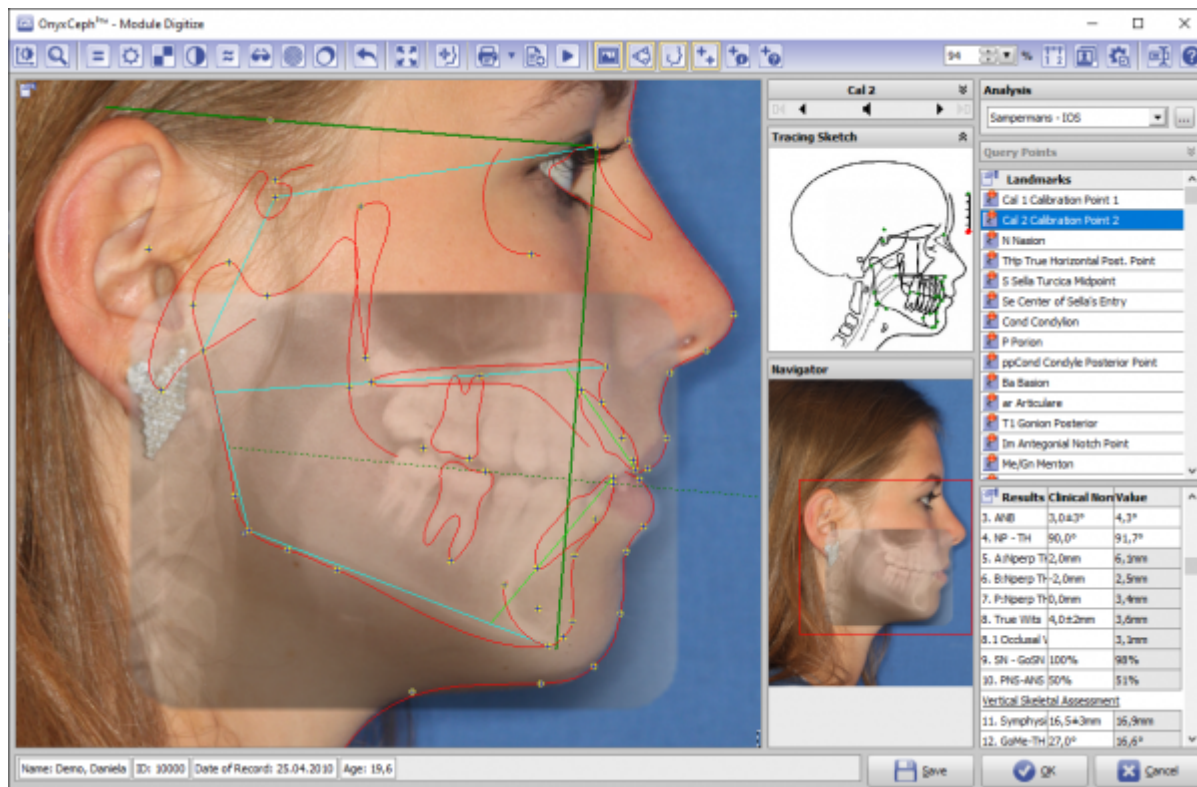
Each analysis comprizes a fixed set of diagnostic variables, which are defined by the position of certain landmarks in the image and, in several cases, from additional input information.

The main function of the module is accordingly to locate the required reference points on the data set and to automatically calculate the analysis results based on them, as well as their graphic representation in the form of tracing and other suitable types of result presentation. Other functions are the output of analysis results or derived information.




First Steps

1. Select analysis
2. Calibrate image
3. Locate landmarks by right click
4. Adjust contour templates
5. Display | print results
6. Save landmarks

Module user interface



Links

	Module Handout Digitize 2D
	2D Analysis Selection
	Calibration of 2D Images

From:

<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

http://onyxwiki.net/doku.php?id=en:digitize_2d&rev=1644219823

Last update: **2022/02/07 08:43**

