

Module Ricketts V.T.O.



Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	✗	✖	✖	✖	✖


✗ not included
 ✖ included by default
 ⓘ optionally included

Module Ricketts V.T.O. [according to the systematics developed by R.M. RICKETTS] can be used to simulate growth and certain orthodontic and/or surgical treatment measures with the help of a lateral ceph tracing in combination with an associated lateral photo in their effect on cephalometric measurements, space discrepancy and soft tissue profile.

First Steps

1. Open module with lateral ceph xray and (optionally) a corresponding lateral photo
2. Define growth period and growth pattern
3. Simulate treatment
4. Analyze cephalometric changes (graphically, numerically)
5. Perform soft tissue simulation: compare before and after
6. Save planning as a project
7. Save modified images (lateral ceph, lateral photo) as findings

Module user interface

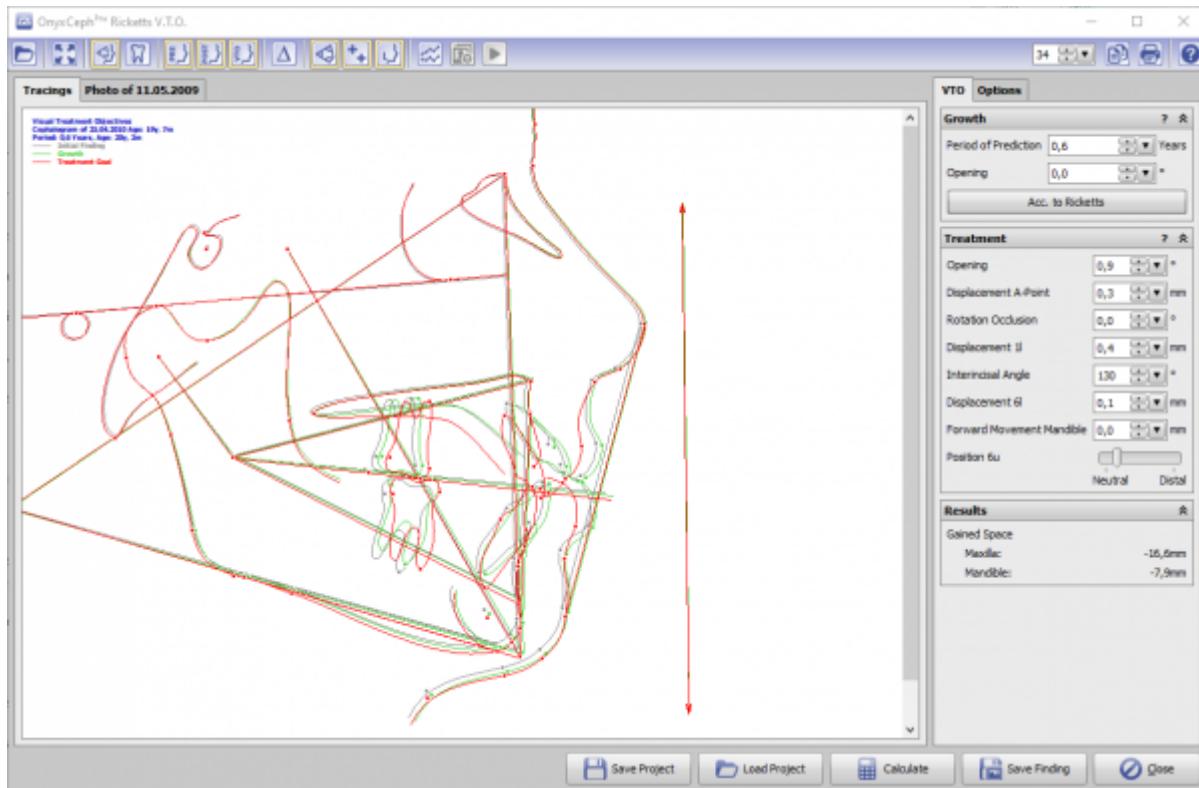

Links



From:
<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
http://onyxwiki.net/doku.php?id=en:ricketts_vto_2d&rev=1676908927

Last update: **2023/02/20 17:02**

