


# Template Editor


## Getting Started

(1) Open empty (new) or existing (selected) print template in print editor,

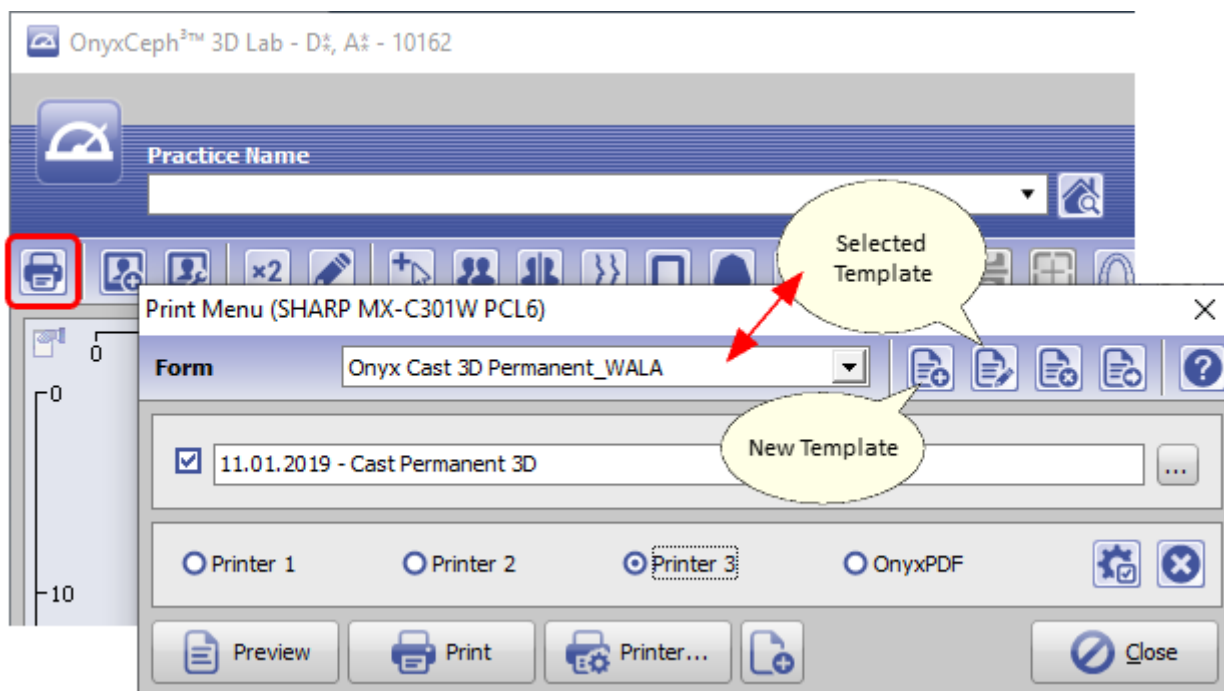

You can do so also via main menu [Options\Systemoptions] on tab [Print]. If you open the editor from the main window directly by icon button [Print], the selected patient and finding is used for the edito-internal previews.

(2) Add or select a print object,

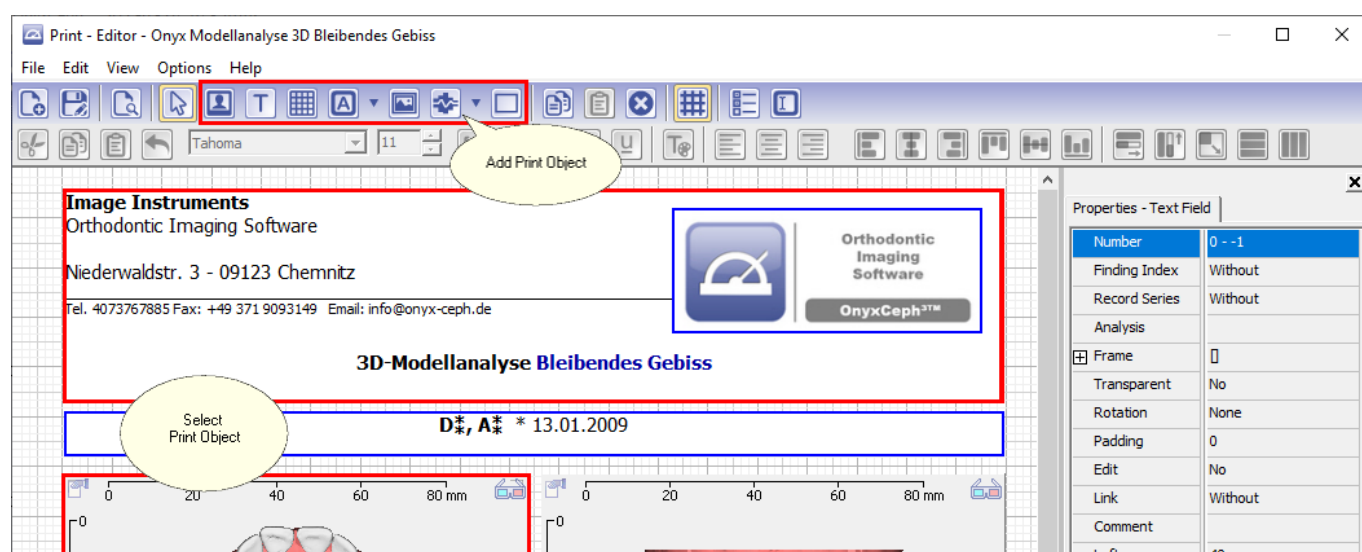

(3) Assign print object position,

... by moving the print object with mouse drag&drop -  
or by changing the numerical position values (top/left) in the Properties dialog

on the right

(4) Print object size

... by moving the print object border with mouse drag&drop -  
or by changing the numerical position values (bottom/right) in the Properties dialog on the right

(5) Print object properties

... by assigning/changing the related property in the Properties dialog on the right

(6) And/or edit objects like text fields and tables by doubleclick.

(7) Finally, save the new/modified template by a unique name before closing the editor.

From:  
<http://onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:  
[http://onyxwiki.net/doku.php?id=en:template\\_editor](http://onyxwiki.net/doku.php?id=en:template_editor)

Last update: **2022/08/23 07:56**

