


# Interface Settings Face Scanner

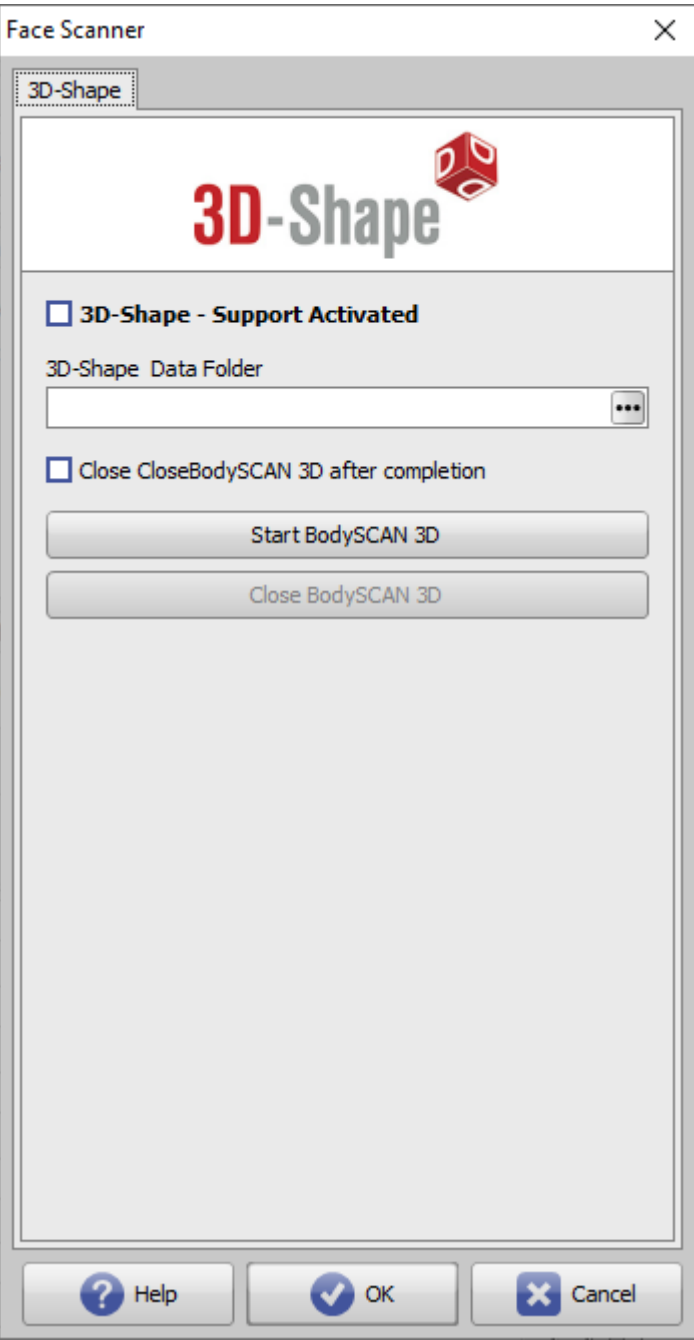

Parameter	Setting
3D Shape - Support Activated	Interface Activation in source panel module <a href="#">Add Image 3D</a>
3D Shape Data Folder	Storage folder for Scan data
Close BodySCAN 3D after completion	Default setting close scan program when scan is completed
Button [Start BodySCAN 3D]	Launch BodySCAN manually - e.g., for testing
Button [Close BodySCAN 3D]	Close BodySCAN manually - e.g., for testing

From:

<http://onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:

<http://onyxwiki.net/doku.php?id=en:3d-shape>

Last update: **2021/03/31 10:10**

