






## Module Cast Adjust




Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability					

Module Cast Adjust allows to re-align the jaw relation of a scan according to different specifications and the attachment of different types of model bases (EU, US, plate, horseshoe bases) after manual definition of a trim line. The study model generated in this way must be saved as a separate finding.

### First Steps

1. Open Module with scan (open shell - two arches in occlusion)
2. Optional: Adjust jaw relation on tab |Articulate|, |Jaw Rotation Axis| or |Align Jaw|
3. On tab |Clip Object| by right clicks, define trim line for upper and lower jaw
4. Clip along trim line Press
5. On tab |Model Base| select type of base tray and merge with scan
6. Save study model as new finding

### Module user interface

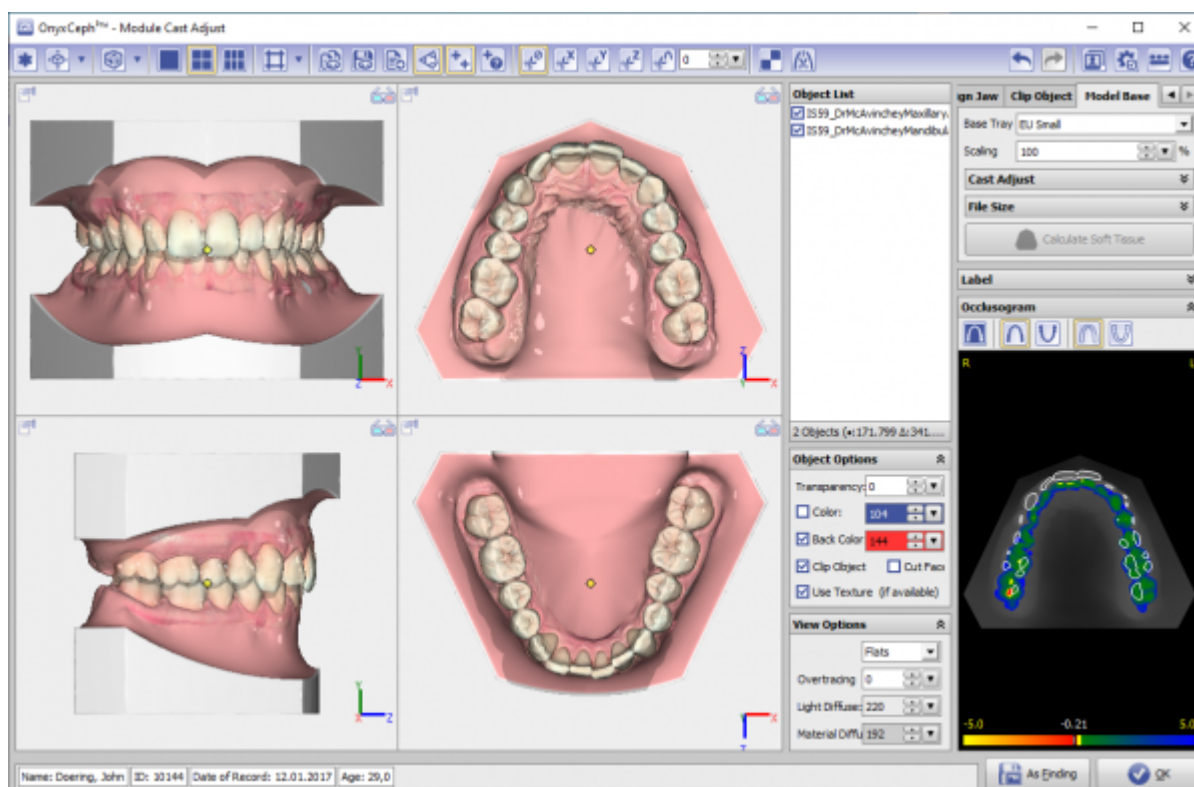

## Links



[https://image-instruments.de/Onyx3/Handout\\_Digitize\\_MA\\_EN.pdf](https://image-instruments.de/Onyx3/Handout_Digitize_MA_EN.pdf)



<https://youtu.be/TAAGkmFIKD4>



<https://drive.google.com/open?id=0B5nAIDRchl5WakjBYXJ2QWNFSXM>

From:

<https://onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:

[https://onyxwiki.net/doku.php?id=en:base\\_tray&rev=1586011014](https://onyxwiki.net/doku.php?id=en:base_tray&rev=1586011014)

Last update: **2020/04/04 16:36**

