

Module Cast Adjust



| Program version | BASIC | 2D PRO | 3D PRO | LAB | OMS |
|-----------------|-------|--------|--------|-----|-----|
| Availability | | | | | |


- nicht
enthalten
- immer
enthalten
- optional
enthalten

Module Cast Adjust allows to re-align the jaw relation of a scan according to different specifications and the attachment of different types of model bases (EU, US, plate, horseshoe bases) after manual definition of a trim line. The study model generated in this way must be saved as a separate finding.

First Steps

1. Open Module with scan (open shell - two arches in occlusion)
2. Optional: Adjust jaw relation on tab |Articulate|,|Jaw Rotation Axis| or |Align Jaw|
3. On tab |Clip Object| by right clicks, define trim line for upper and lower jaw
4. Clip along trim line Press
5. On tab |Model Base| select type of base tray and merge with scan
6. Save study model as new finding

Module user interface

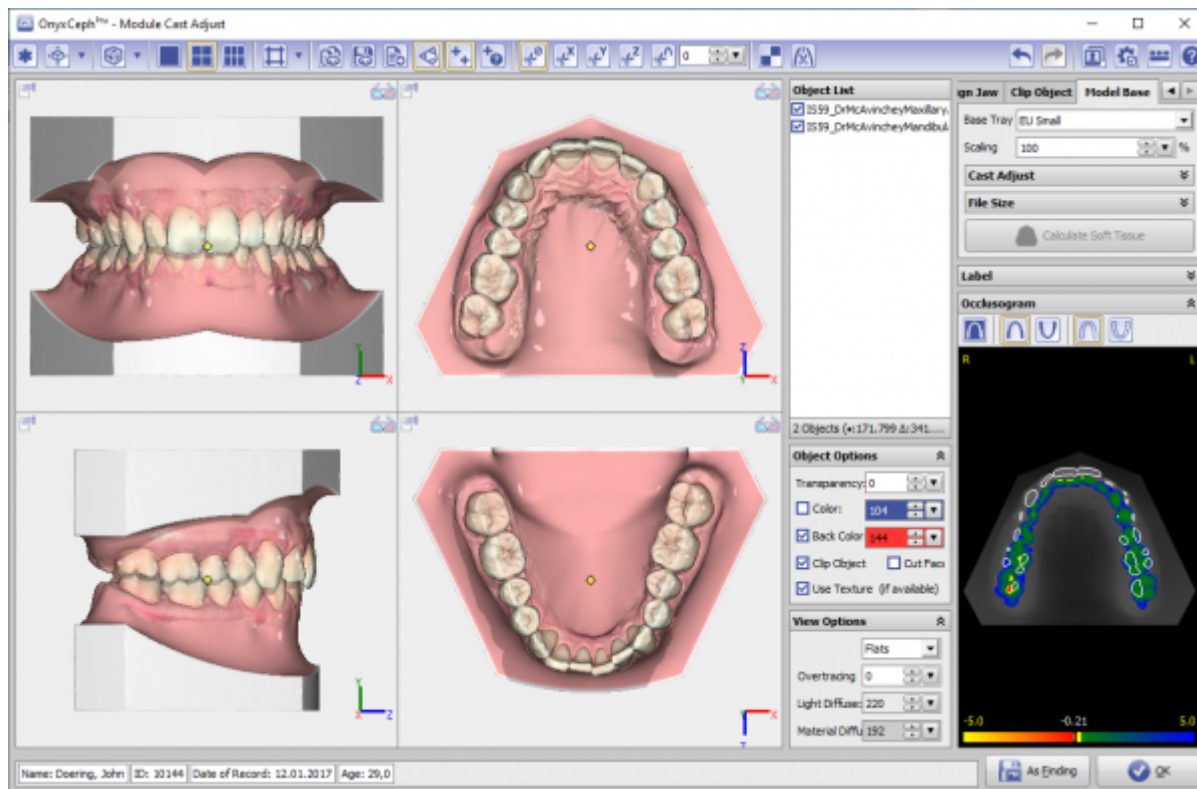

Links



Module Handout Cast Adjust



<https://youtu.be/TAAGkmFIKD4>



<https://drive.google.com/open?id=0B5nAIDRchI5WakJBYXJ2QWNFSXM>

From:

<https://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

https://onyxwiki.net/doku.php?id=en:base_tray&rev=1644236612

Last update: **2022/02/07 13:23**

