

# Module Combine 2D




Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability					

Module Combine 2D can combine two images with different but overlapping content. The second image behind the first image is shown completely or in a self-defined region of interest using common reference points with selectable transparency. In addition, the tracing of one of the two images can be overlaid.

## First Steps

1. Open module with 2 images
2. Select or digitize correspondin landmarks
3. Define region of interest
4. Switch to tab |Combine|
5. Optionally, modify superimposition
6. Save result as project and/or as finding

## Module user interface

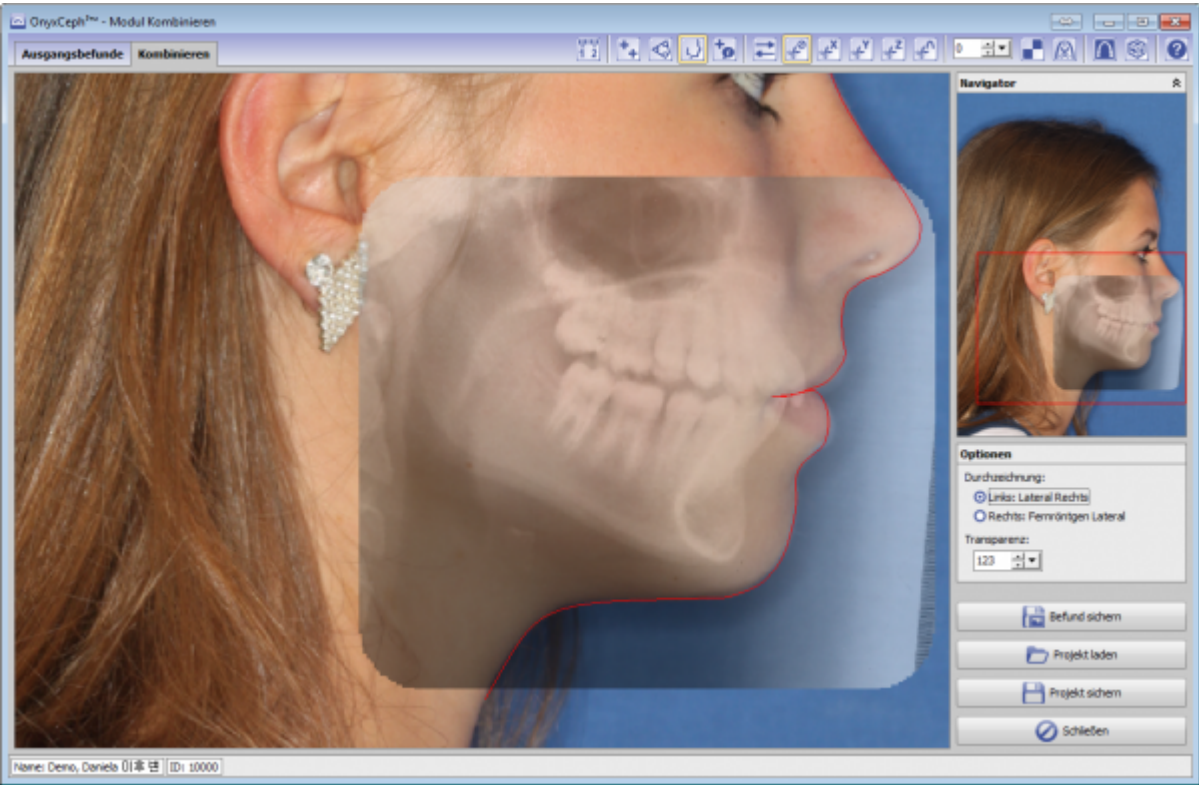

## Links



From:

<https://onyxwiki.net/> - [OnyxCeph<sup>3</sup>™ Wiki]

Permanent link:

[https://onyxwiki.net/doku.php?id=en:combine\\_2d&rev=1586009762](https://onyxwiki.net/doku.php?id=en:combine_2d&rev=1586009762)

Last update: **2020/04/04 16:16**

