

Module TADmatch™



Program version	BASIC	2D PRO	3D PRO	LAB	OMS
Availability	✗	✗	📦	📦	📦

- ✗

enthalten

nicht
- 📦

enthalten

immer
- 📦

enthalten


optional

The module TADmatch allows a virtual placement of anchorage devices from a library included in OnyxCeph with the intention of calculating models which are either used for transferring the placement to the patient or manufacturing the appliance on the model.

First Steps

1. Combine relevant findings in the module Register3D, save new finding
2. In new finding, select item from library and place on the model
3. Compute positioning or manufacturing model
4. Export models using the object list
5. Save as project and new finding

Module user interface

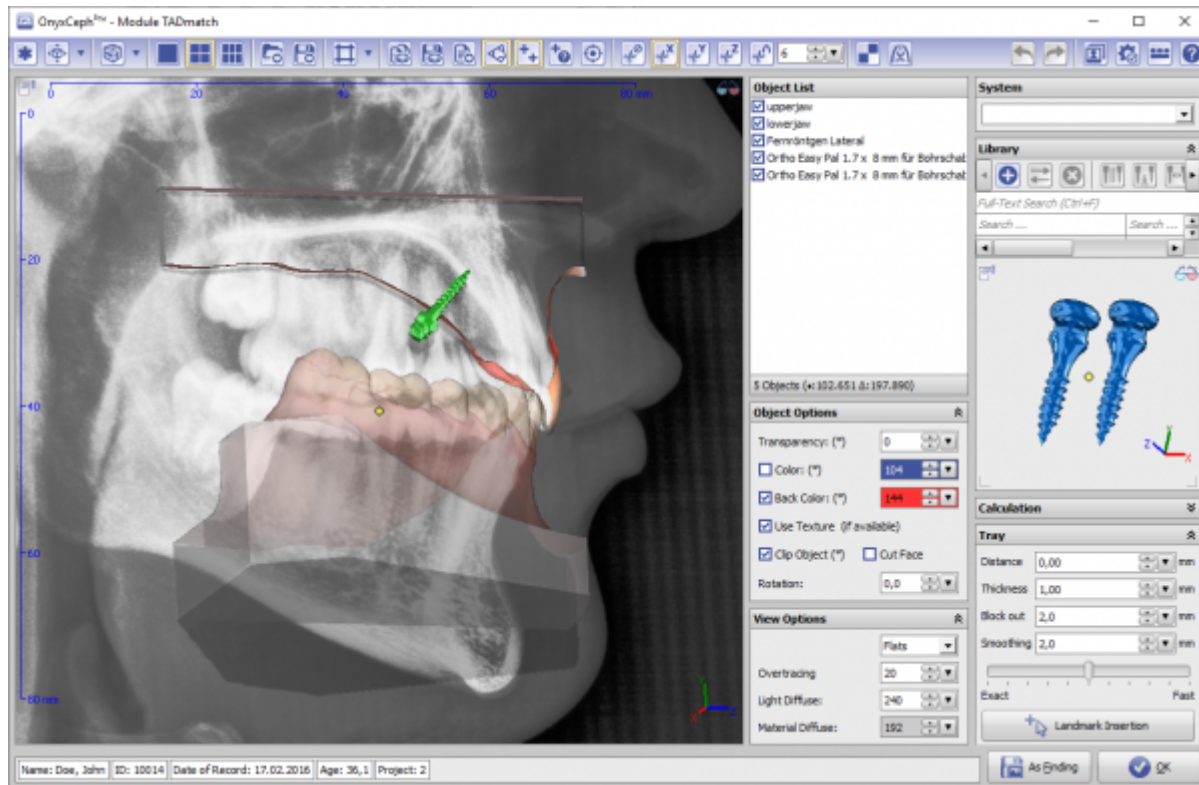

Links



From:
<https://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
https://onyxwiki.net/doku.php?id=en:tadmatch_3d&rev=1644236968

Last update: **2022/02/07 13:29**

