




This page is not fully translated, yet. Please help completing the translation.
(remove this paragraph once the translation is finished)

[Medienimport]

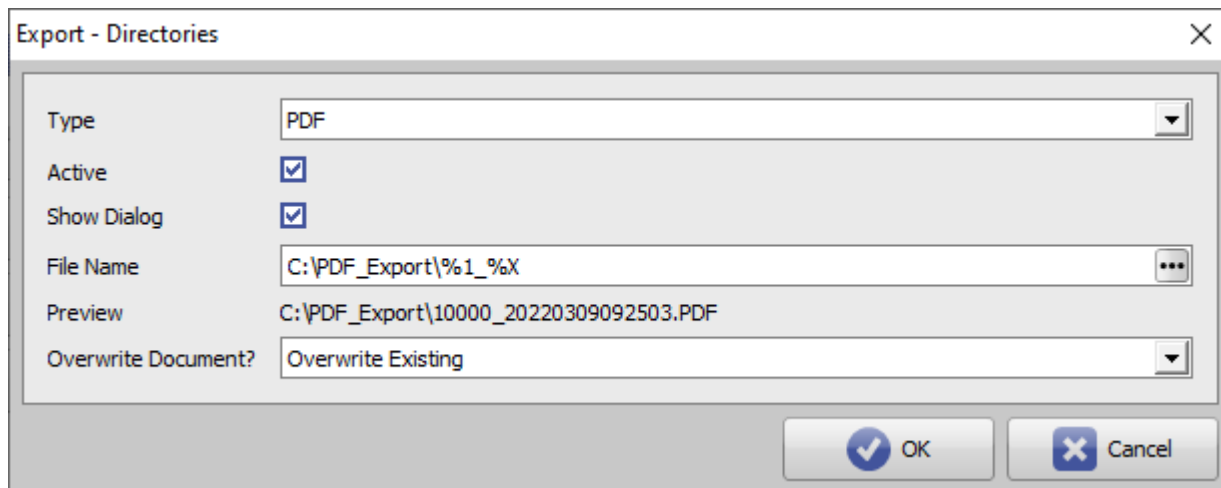

The image shows a Windows-style dialog box titled "Export - Directories". It contains several settings for exporting files:

- Type:** A dropdown menu set to "PDF".
- Active:** A checked checkbox.
- Show Dialog:** A checked checkbox.
- File Name:** A text field containing "C:\PDF_Export\%1_%X" with a browse button (three dots) to its right.
- Preview:** A text field showing "C:\PDF_Export\10000_20220309092503.PDF".
- Overwrite Document?:** A dropdown menu set to "Overwrite Existing".

At the bottom right, there are two buttons: "OK" (with a blue checkmark icon) and "Cancel" (with a blue X icon).

Edit: Image Type

General

Image Type Name: Lateral Ceph

Image Type Category: X-Ray

Token: XRCR

Thumbnail: 

Load...

Delete

Image Acquisition

Maximum Image Size: 4096

Image Aspect Ratio (%): 64,84

Horizontal Guideline (%): 40

Vertical Guideline (%): 80

Rotate Image: 0°

Fit Image to Frame: Fit All

Standard...

☐ Mirror Image Horizontally

☐ Mirror Image Vertically

Tracing

Tracing Type: Ceph Lateral

☐ Trace as Flipped

Analyzing Methods



Slide View

View: Riyahd

Printing Form: KSA MOH Ceph

Help

OK

Cancel

From:
<http://www.onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
<http://www.onyxwiki.net/doku.php?id=en:mediaimport>

Last update: **2022/03/09 10:13**

