

Scan Import Medit


As of release OnyxCeph³™ 3.2.157 (438), direct scan data import from Medit Link Online Case Portal is supported. Here, there are 2 ways to transfer Medit scans to OnyxCeph³™ :

POST	Transfer a scan from Medit Link via the OnyxCeph ³ ™ Integration App to OnyxCeph ³ ™ Module Add 3D Image
GET	Calling Medit Link from OnyxCeph ³ ™ module Add 3D image and retrieving captured scan datasets.

While any scans can be transferred to OnyxCeph³™ using the **POST** method, only the scan data created by calling OnyxCeph³™ can be transferred using the **GET** method. Prerequisite for the use of both methods is the installation of the OnyxCeph³™ Integration App in the Medit Link Online Case Portal.


Settings Medit Link

Install the OnyxCeph³™ Integration App from the Medit Link App Store:

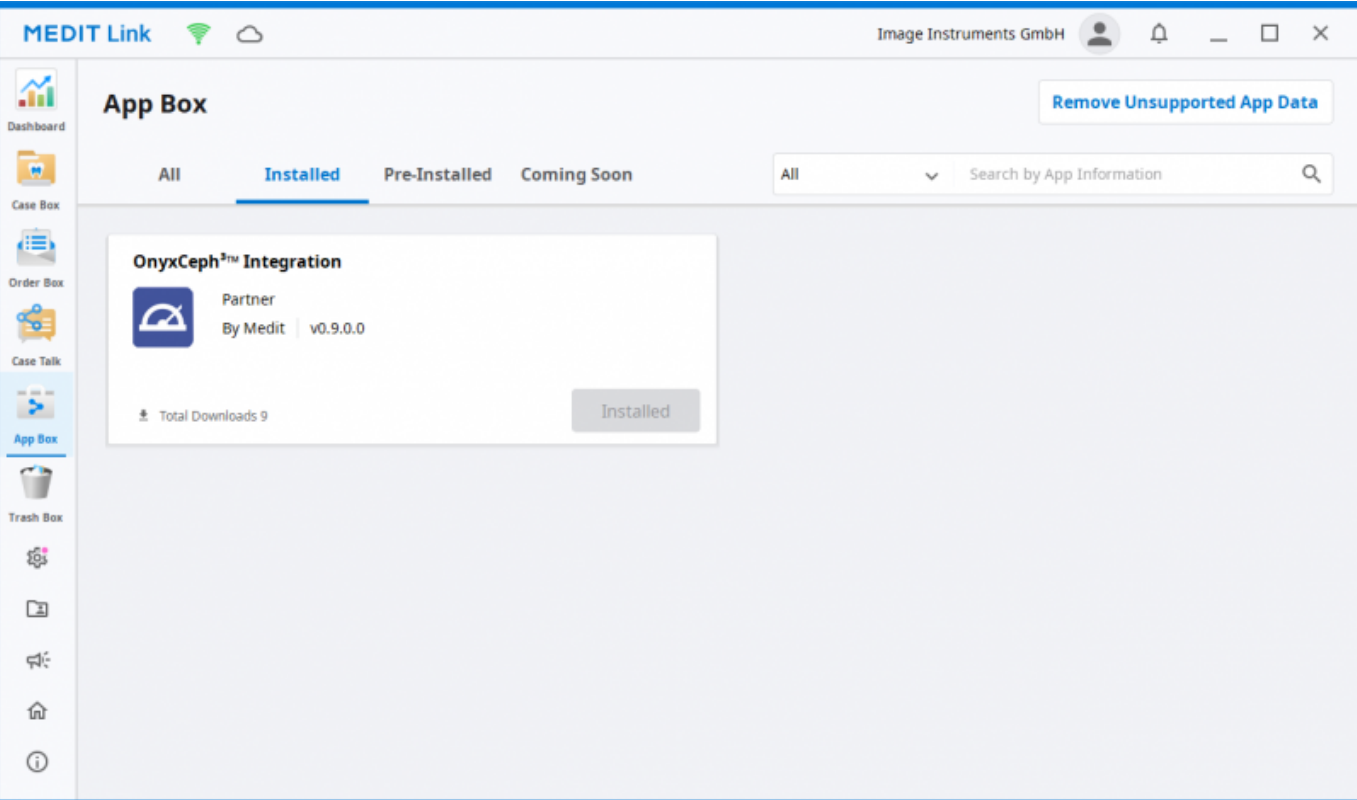


Settings OnyxCeph³™

1. Open main menu Options[System Options on [tab |interfaces|](#)].
2. In the [Import](#) panel, click on [\[Model - i/o Scanner\]](#) button.
3. Open [Medit Link](#) tab
4. Select program folder medit_Link.exe
5. Enter offset values for rotation correction (if necessary)
6. Activate checkbox [Medit Link - support active]

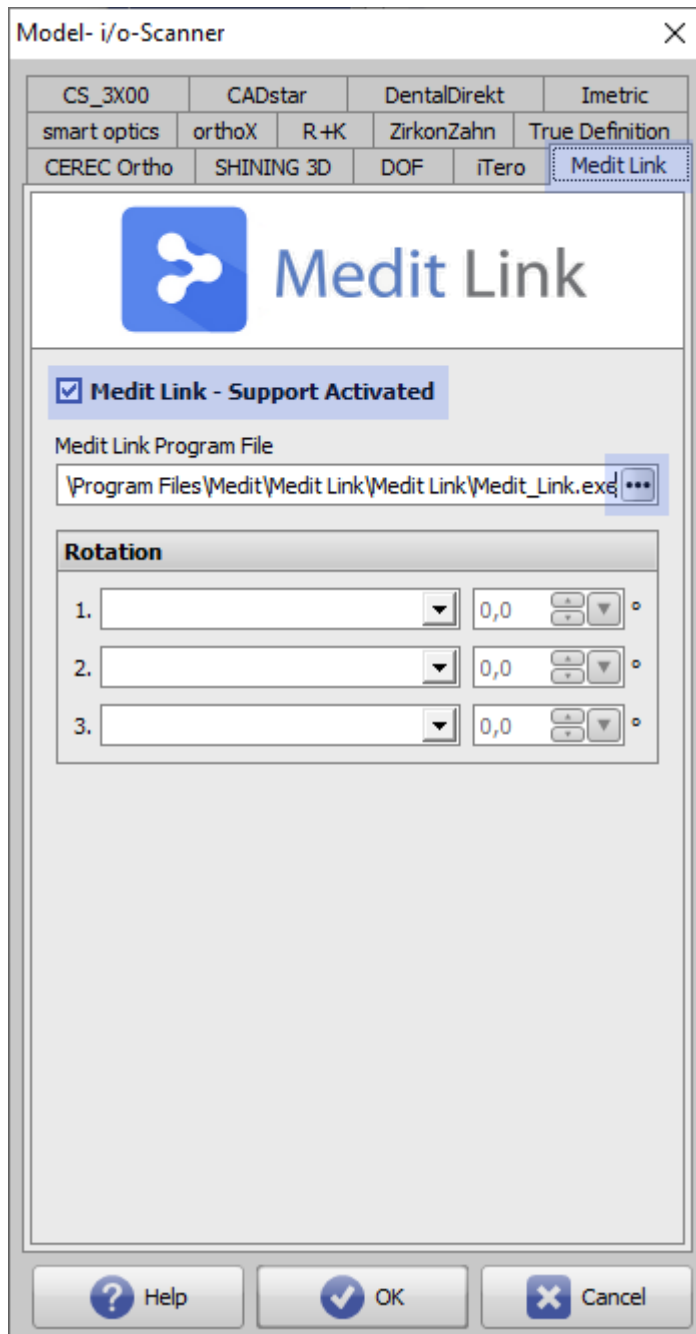

Model- i/o-Scanner

CS_3X00	CADstar	DentalDirekt	Imetric
smart optics	orthoX	R+K	ZirkonZahn
CEREC Ortho	SHINING 3D	DOF	iTero
Medit Link			

 **Medit Link**

☒ **Medit Link - Support Activated**


Medit Link Program File
\\Program Files\\Medit\\Medit Link\\Medit Link\\Medit_Link.exe

Rotation


1.	<input type="text"/>	<input type="text" value="0,0"/>	<input type="button" value="+"/>	<input type="button" value="-"/>	<input type="button" value="°"/>
2.	<input type="text"/>	<input type="text" value="0,0"/>	<input type="button" value="+"/>	<input type="button" value="-"/>	<input type="button" value="°"/>
3.	<input type="text"/>	<input type="text" value="0,0"/>	<input type="button" value="+"/>	<input type="button" value="-"/>	<input type="button" value="°"/>

Transfer Scan (POST)


1. Open Medit Link on local OnyxCeph^{3™} client computer
2. Select patient | case | scan
3. Click icon button [Export to OnyxCeph]

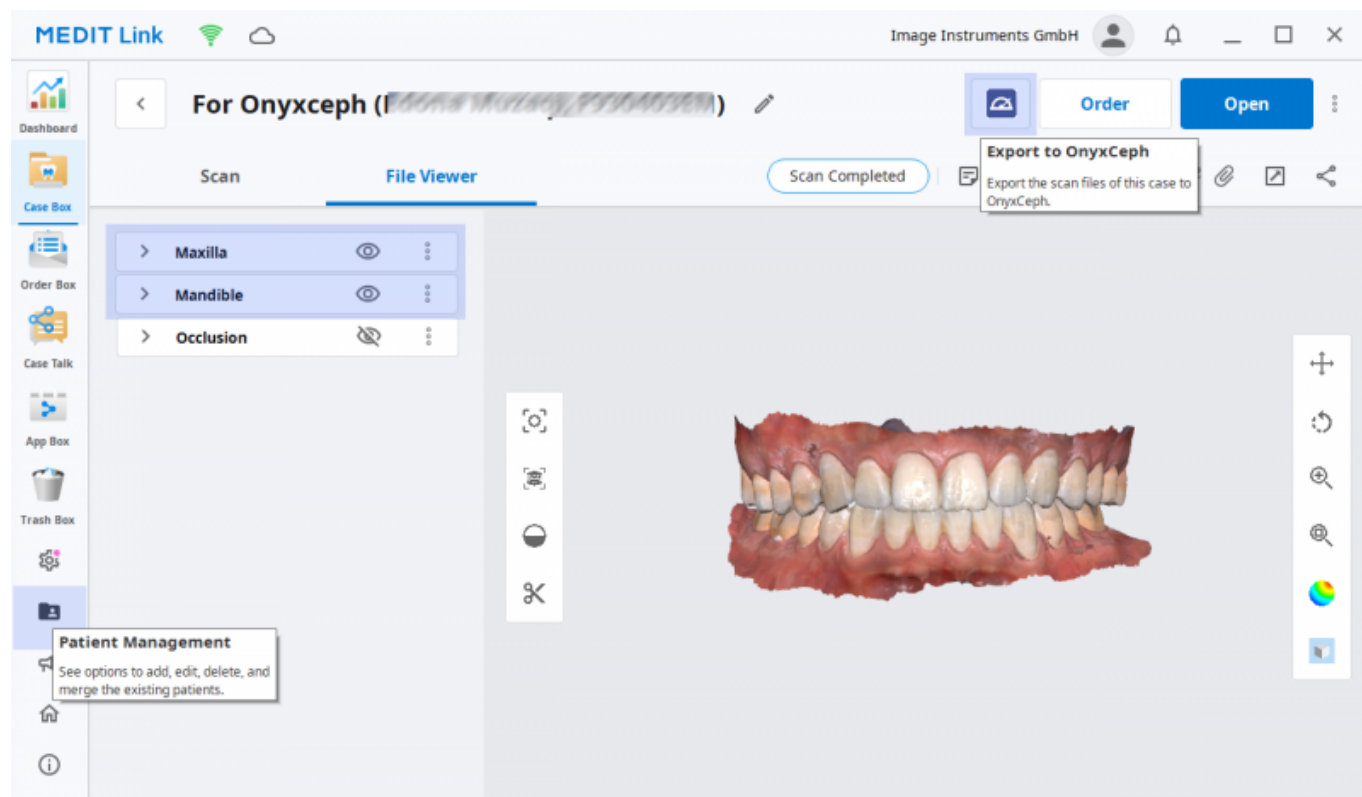

4. If the scan is not associated with the current login, a warning will appear:

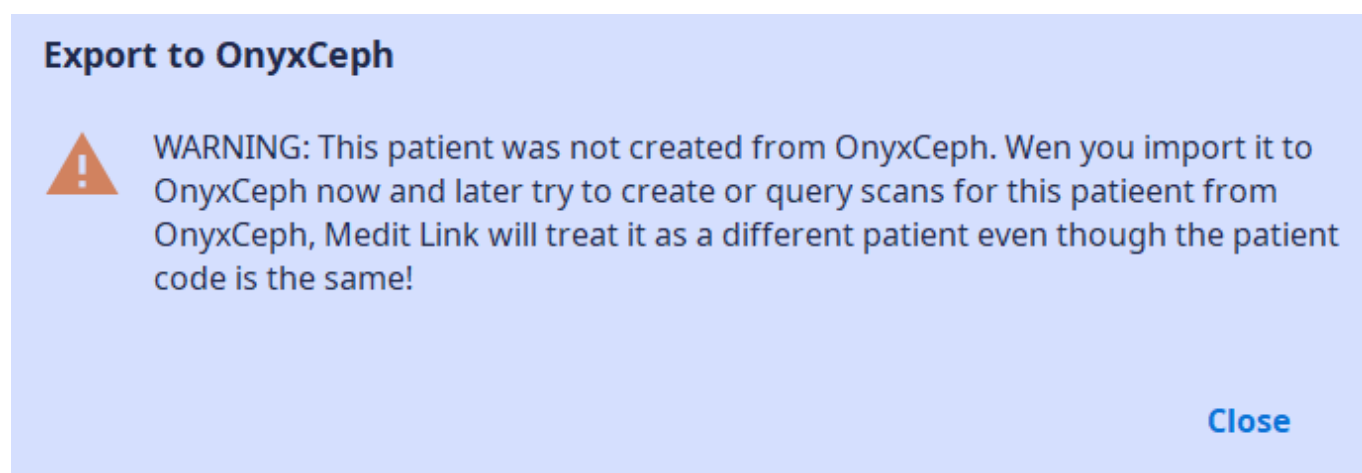


5. If this warning is closed, the scan data is transferred to module [Add image 3D](#) and can be classified there, aligned to the patient and further processed according to the known workflow.

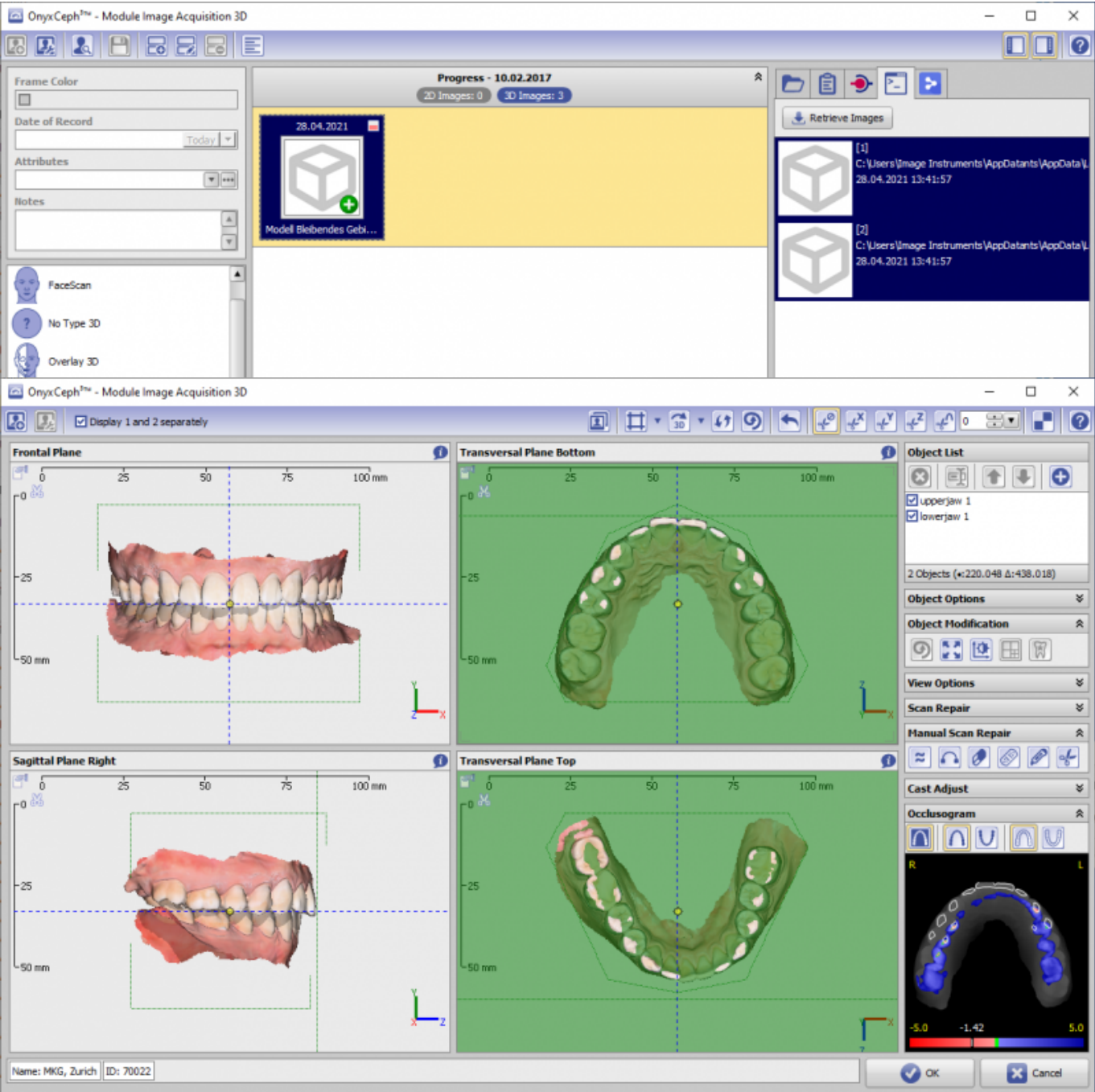


Request Scan (GET)


1. Select patient record in OnyxCeph³™
2. Open module [Bild hinzufügen 3D](#)
3. In the source panel, select tab [Medit Link]
4. Click button [Image Acquisition...] klicken

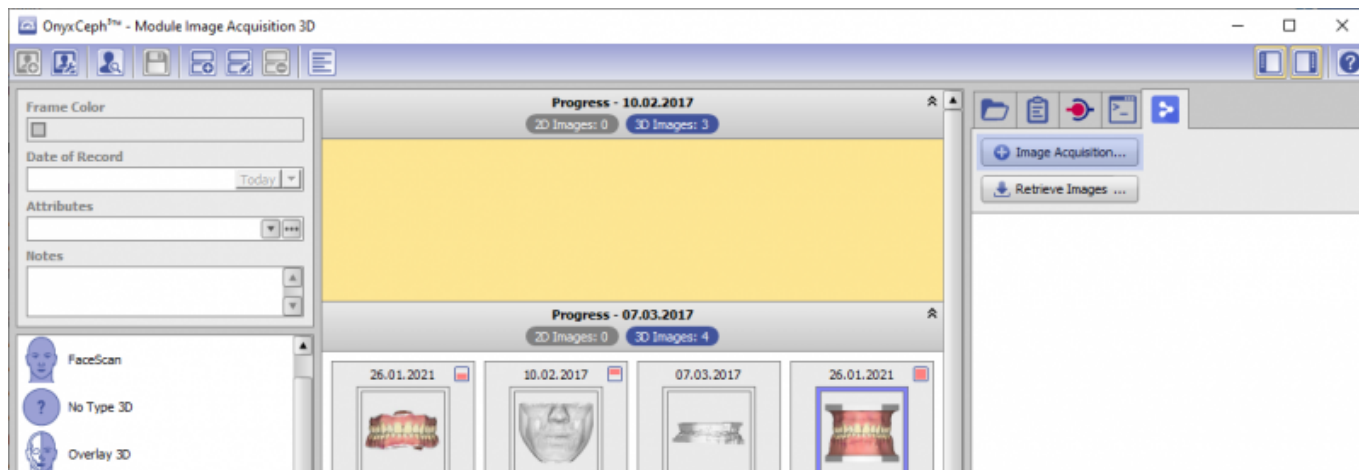


5. In Medit Link, acquire or select scan data

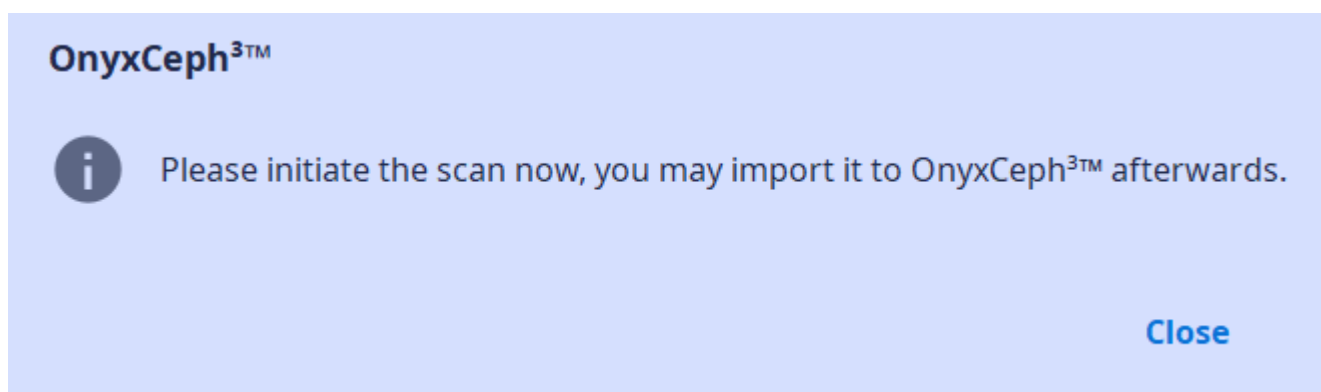

6. List by button [Retrieve Images...] and import, classify, align and further process the scan data according to the known workflow.

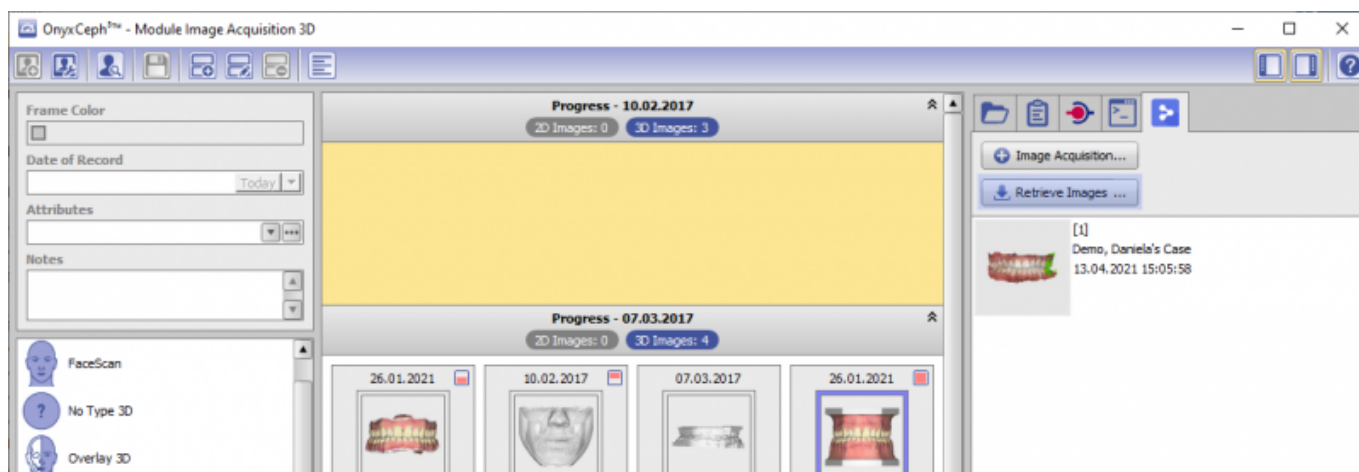

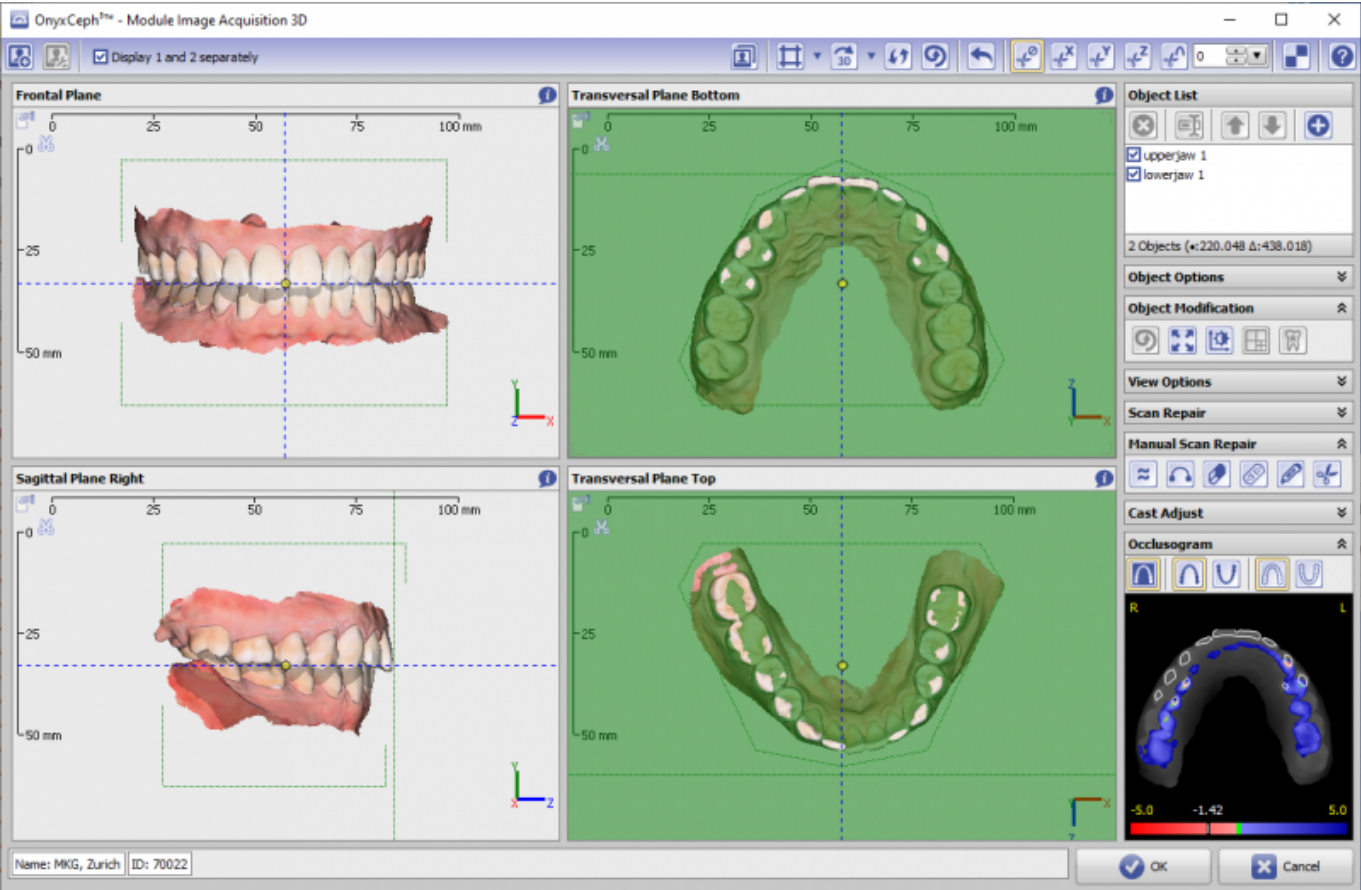

From:

<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

<http://onyxwiki.net/doku.php?id=en:meditimport&rev=1620995761>

Last update: **2021/05/14 14:36**

