

Scan Import Medit


As of release OnyxCeph^{3™} 3.2.157 (438), direct scan data import from Medit Link Online Case Portal is supported. Here, there are 2 ways to transfer Medit scans to OnyxCeph^{3™}:

POST	Medit Link ► OnyxCeph ^{3™}
GET	OnyxCeph ^{3™} ◀ Medit Link

While any scans can be transferred to OnyxCeph^{3™} using the **POST** method, only the scan data assigned to OnyxCeph^{3™} can be transferred using the **GET** method. Prerequisite for the use of both methods is the installation of the OnyxCeph^{3™} Integration App in the Medit Link Online Case Portal.


Settings Medit Link

Install the OnyxCeph³™ Integration App from the Medit Link App Store:

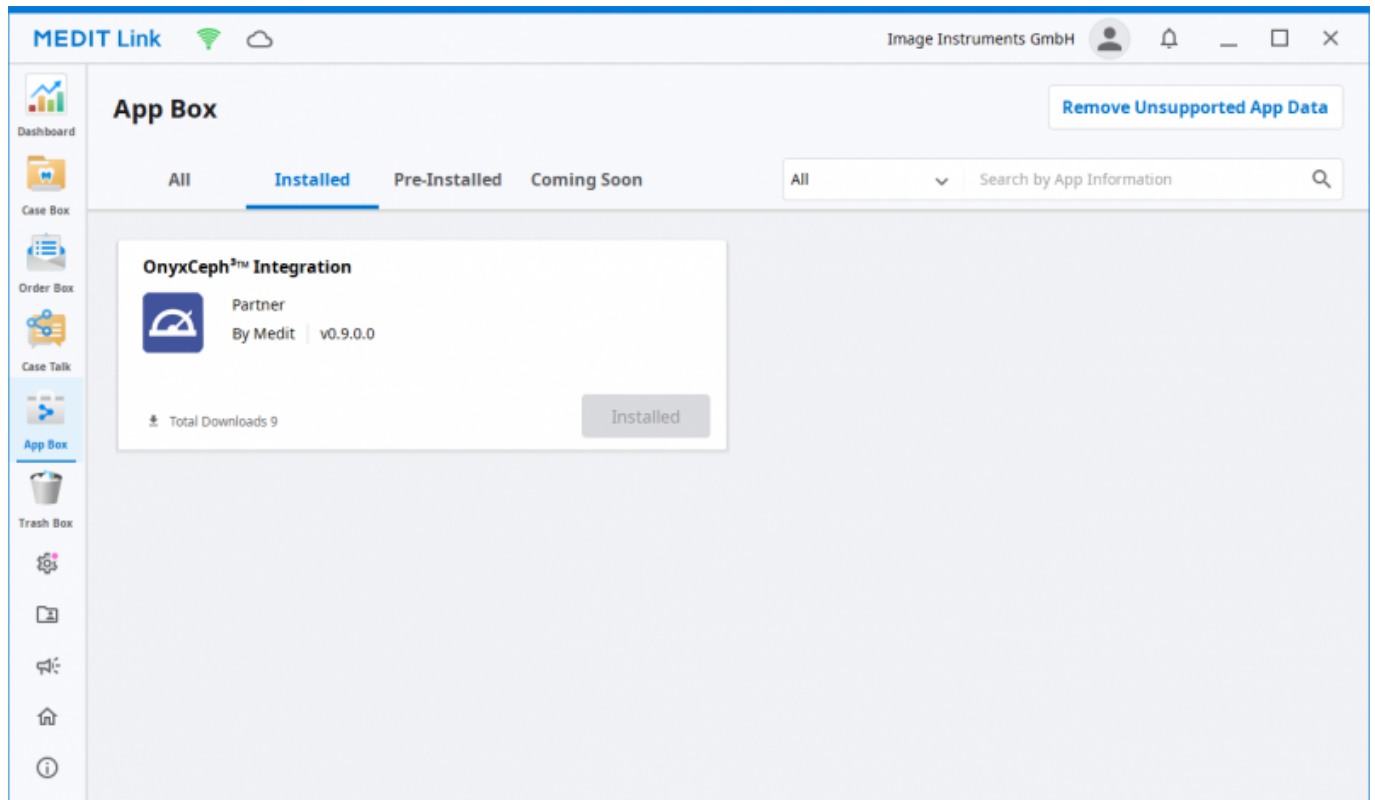


Settings OnyxCeph³™

1. Open main menu Options[System Options on [tab |interfaces|](#)\\.
2. In the [Import](#) panel, click on [\[Model - i/o Scanner\]](#) button.
3. Open [Medit Link|](#) tab
4. Select program folder medit_Link.exe
5. Enter offset values for rotation correction (if necessary)
6. Activate checkbox [Medit Link - support active]

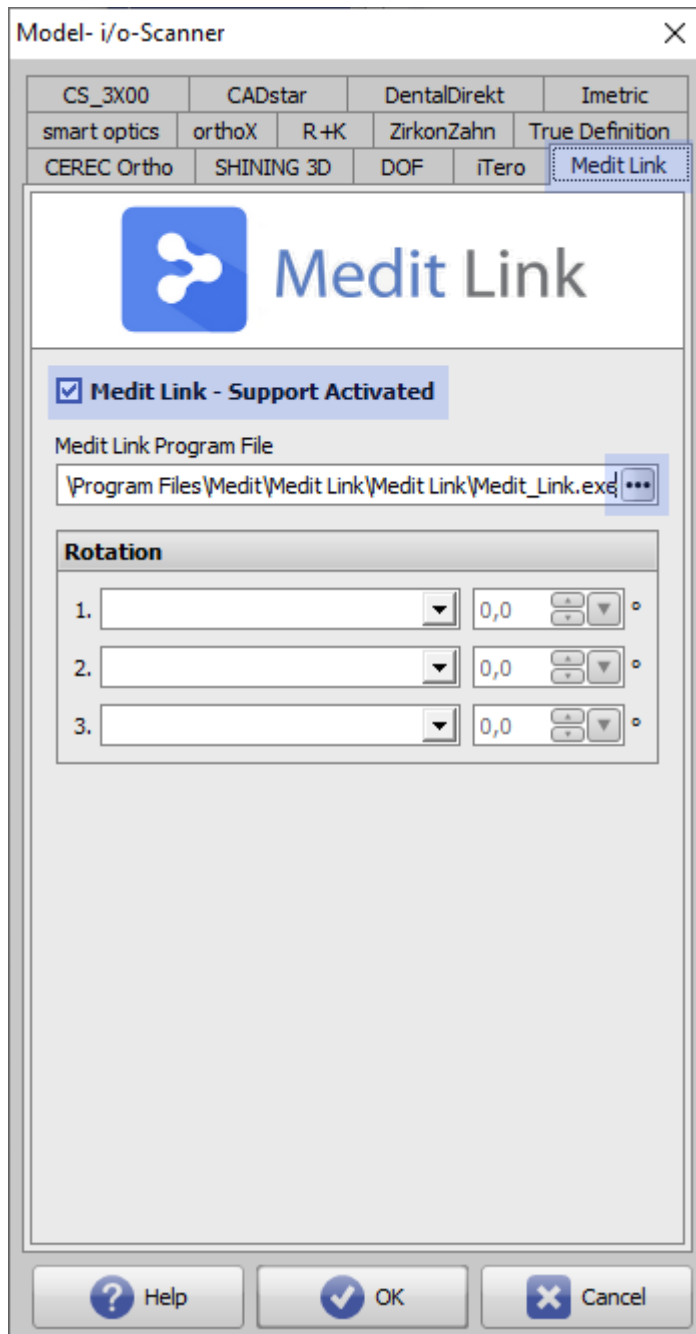


Transfer Scan (POST)

1. Open Medit Link on local OnyxCeph^{3™} client computer
2. Select patient | case | scan ①
3. Click icon button [Export to OnyxCeph]


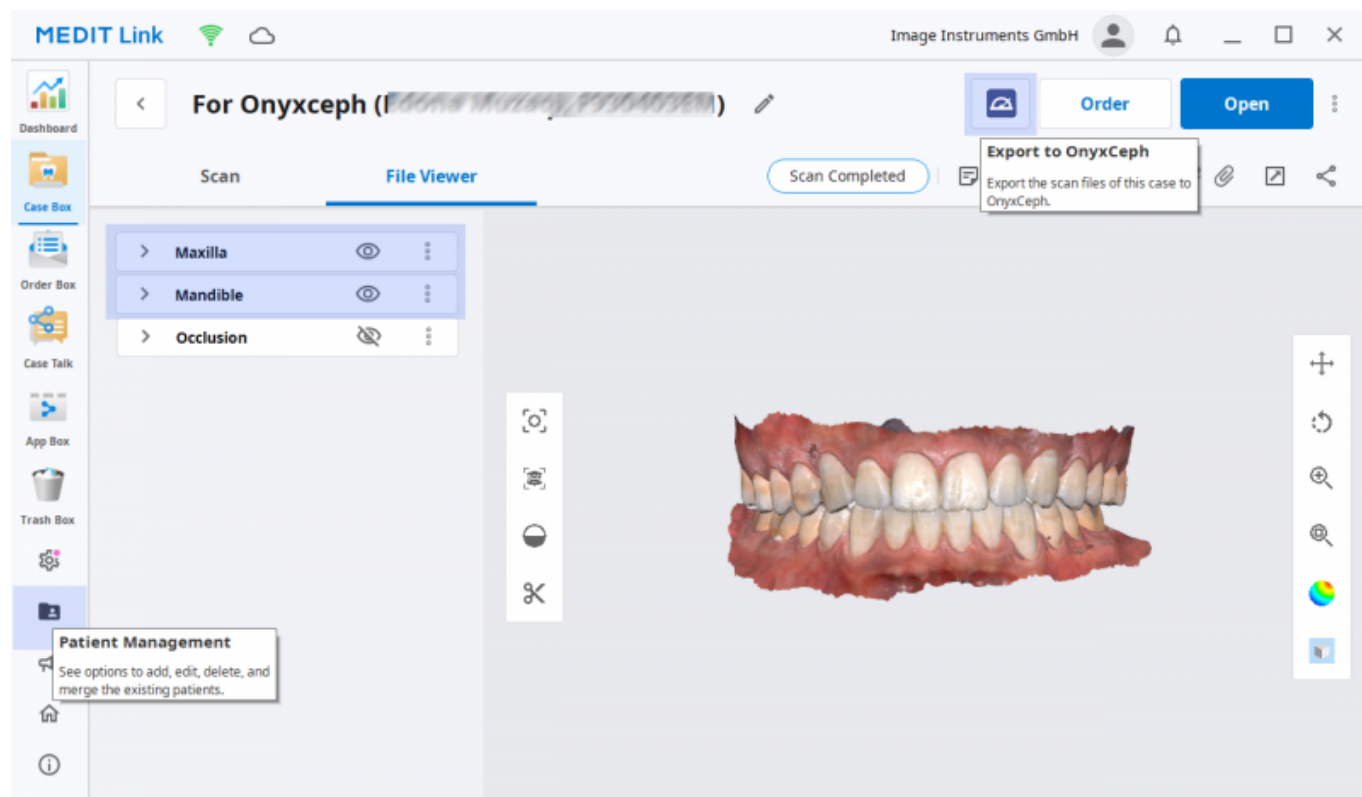


① Patient muss in der Medit Link Patientenverwaltung mit allen Stammdaten angelegt sein!

4. If the scan is not associated with the current login, a warning will appear:

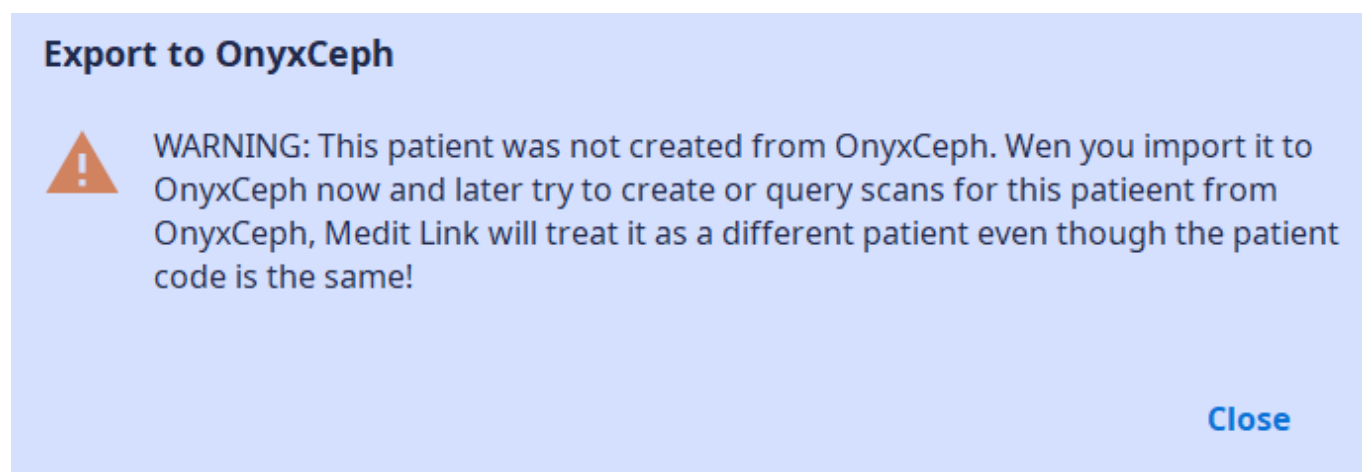


5. If this warning is closed, the scan data is transferred to module [Add image 3D](#) and can be classified there, aligned to the patient and further processed according to the known workflow.

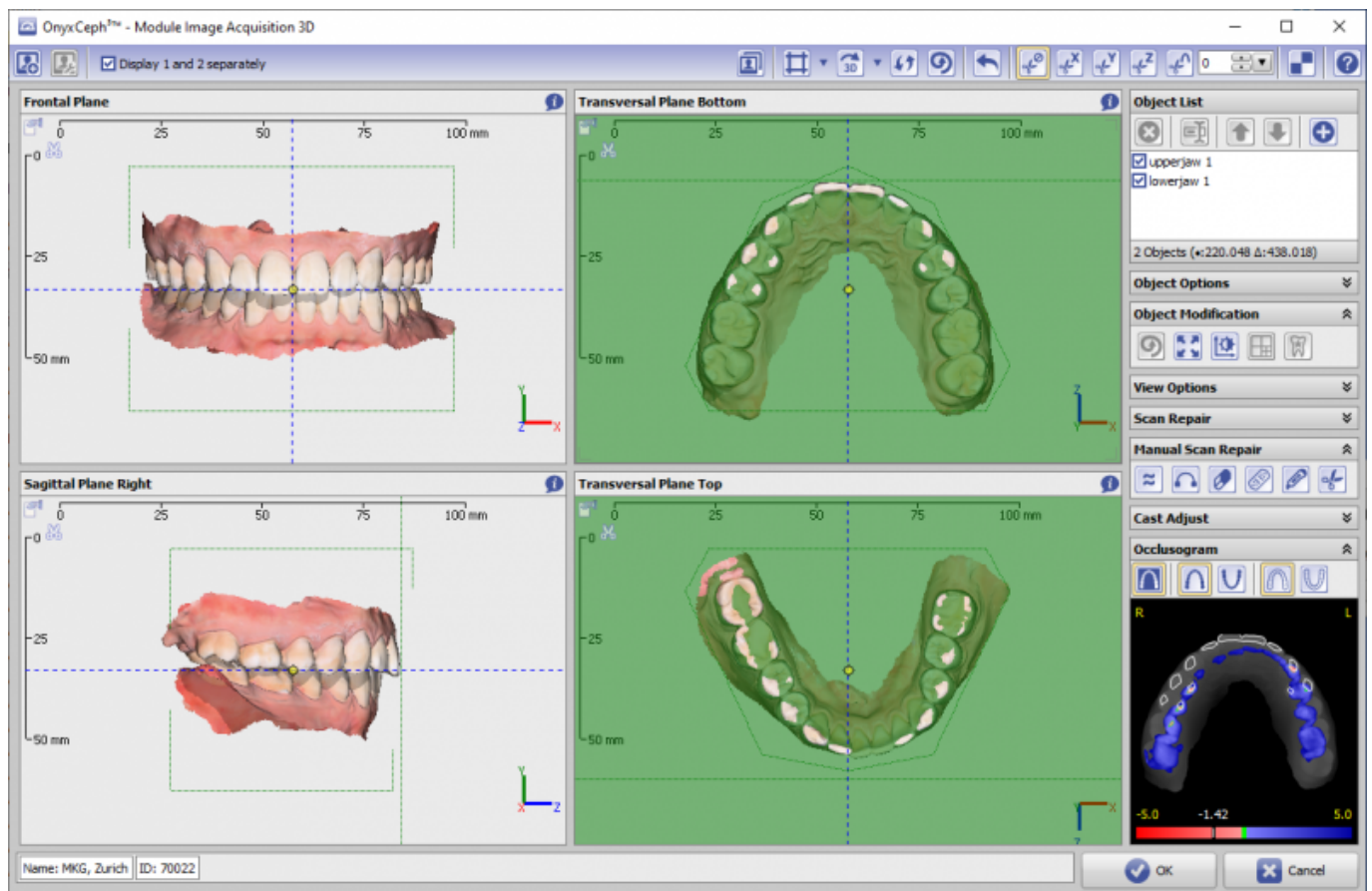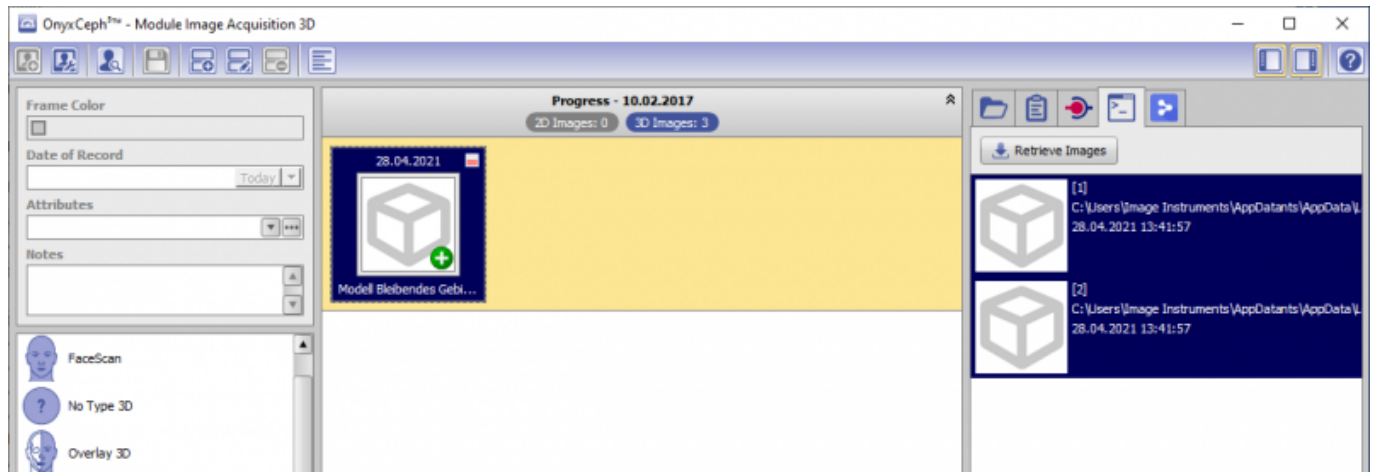


Request Scan (GET)


1. Select patient record in OnyxCeph³™
2. Open module [Bild hinzufügen 3D](#)
3. In the source panel, select tab [Medit Link]
4. Click button [Image Acquisition...] klicken

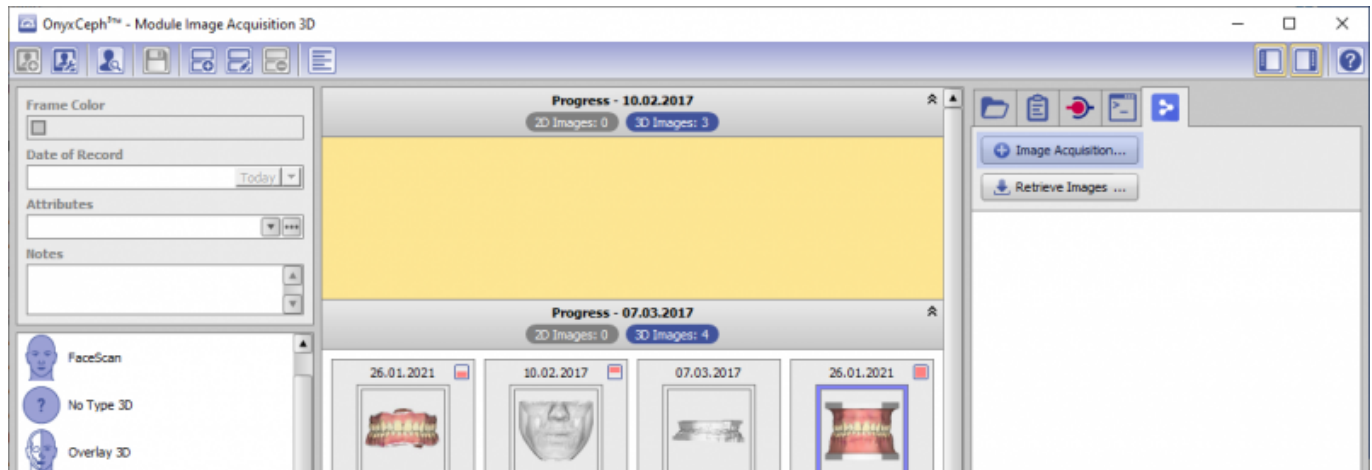


5. In Medit Link, acquire or select scan data

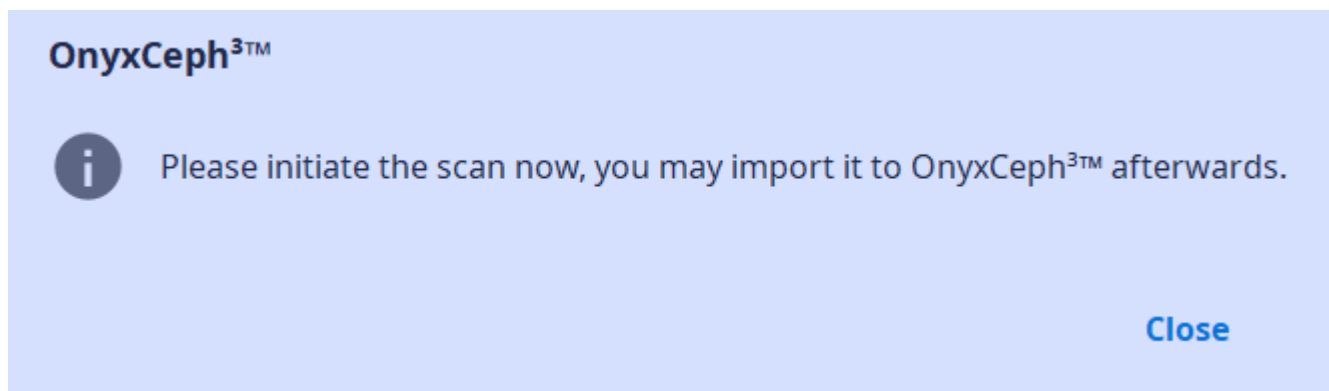

6. List by button [Retrieve Images...] and import, classify, align and further process the scan data according to the known workflow.

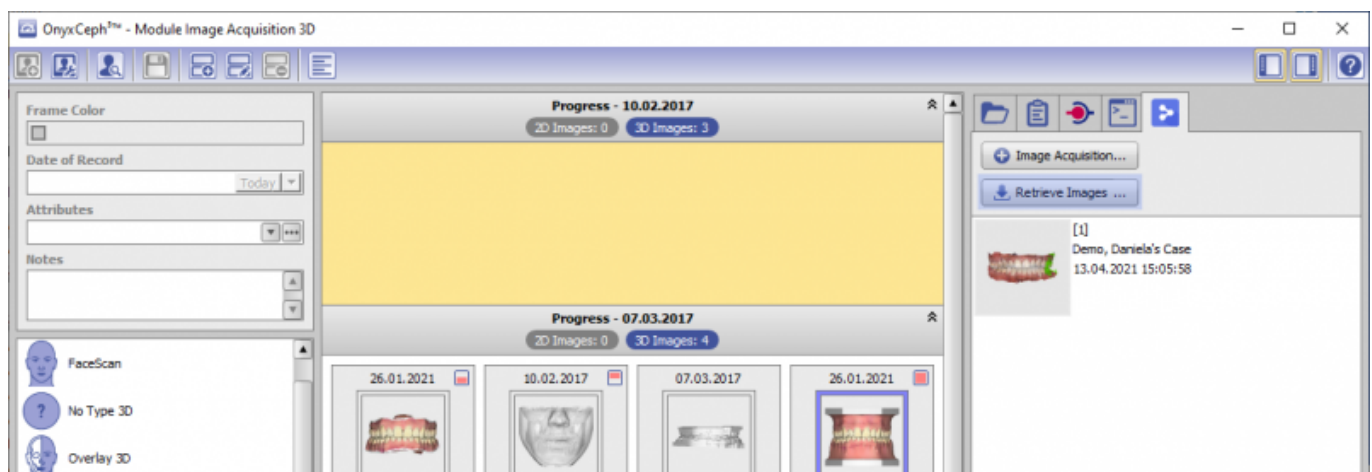

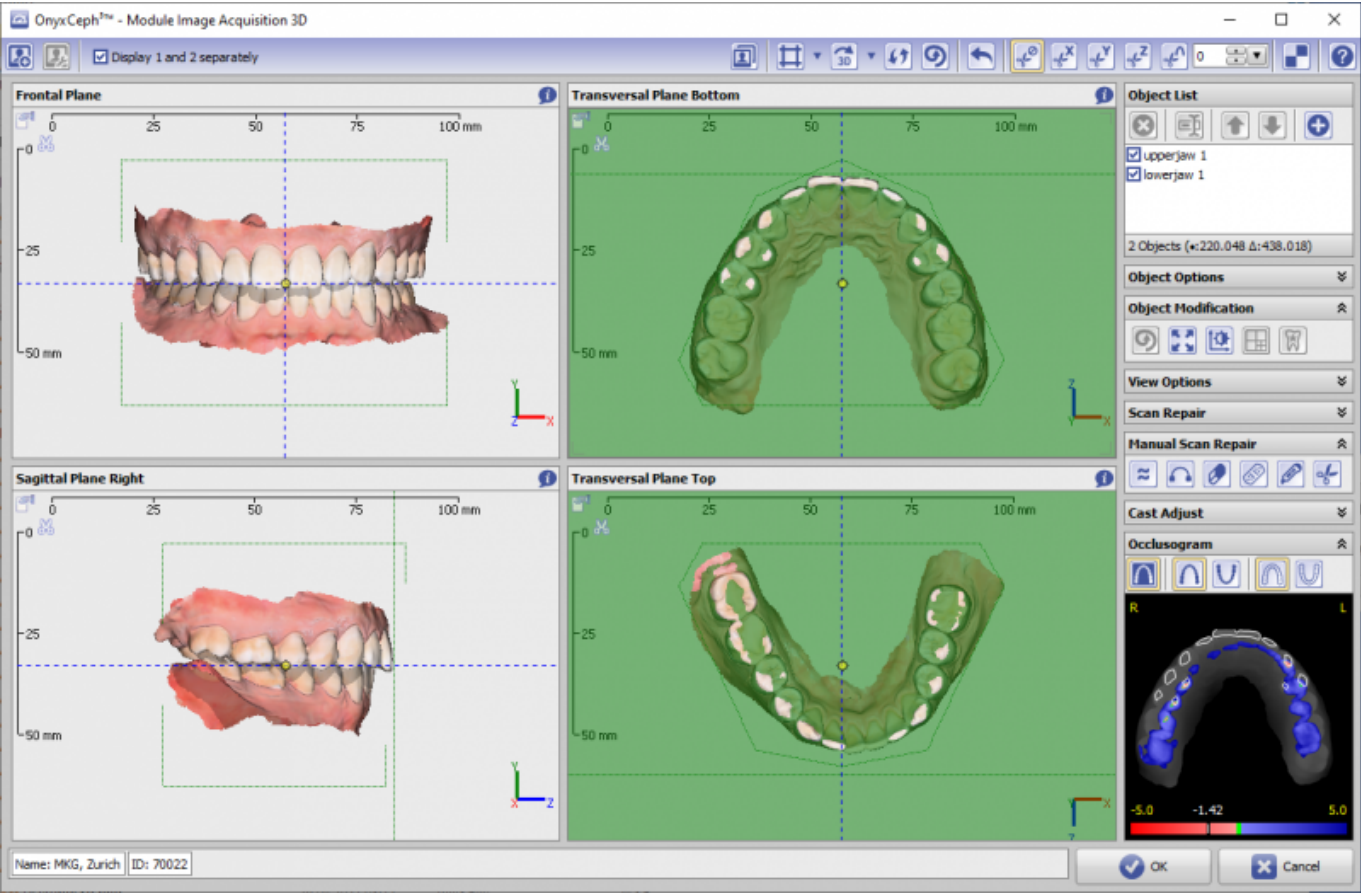

From:

<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

<http://onyxwiki.net/doku.php?id=en:meditimport&rev=1621521949>

Last update: **2021/05/20 16:45**

