

Scan Import Medit


As of release OnyxCeph³™ 3.2.157 (438), direct scan data import from Medit Link Online Case Portal is supported. Here, there are 2 ways to transfer Medit scans to OnyxCeph³™ :

POST	Medit Link ► OnyxCeph ³ ™
GET	OnyxCeph ³ ™ ◄ Medit Link

While any scans can be transferred to OnyxCeph³™ using the **POST** method, only the scan data assigned to OnyxCeph³™ can be transferred using the **GET** method. Prerequisite for the use of both methods is the installation of the OnyxCeph³™ Integration App in the Medit Link Online Case Portal.


Settings Medit Link

Install the OnyxCeph³™ Integration App from the Medit Link App Store. MeditLink version required for this is 2.5.0:

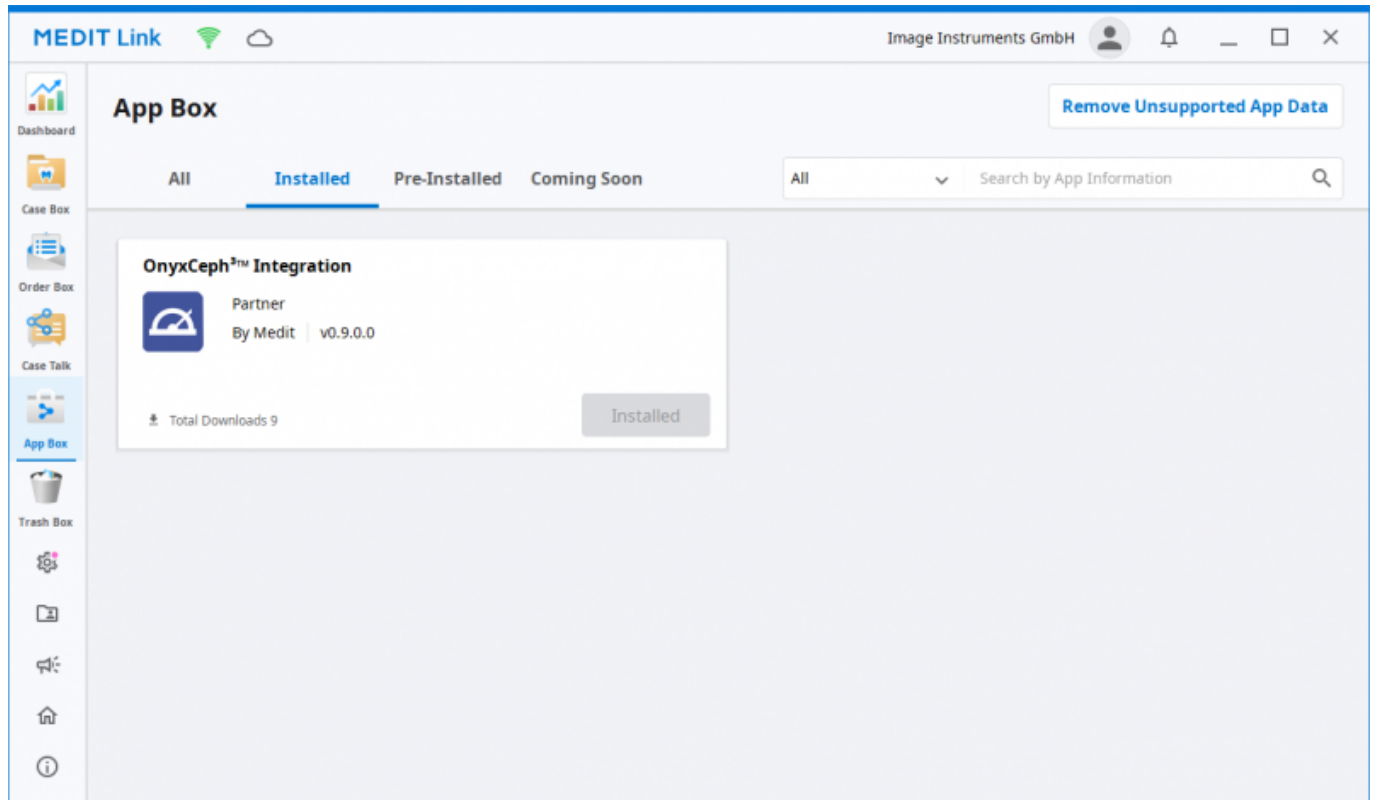


Settings OnyxCeph³™

1. Open main menu Options[System Options on [tab |interfaces|](#)\\.
2. In the [Import](#) panel, click on [\[Model - i/o Scanner\]](#) button.
3. Open [Medit Link|](#) tab
4. Select program folder medit_Link.exe
5. Enter offset values for rotation correction (if necessary)
6. Activate checkbox [Medit Link - support active]

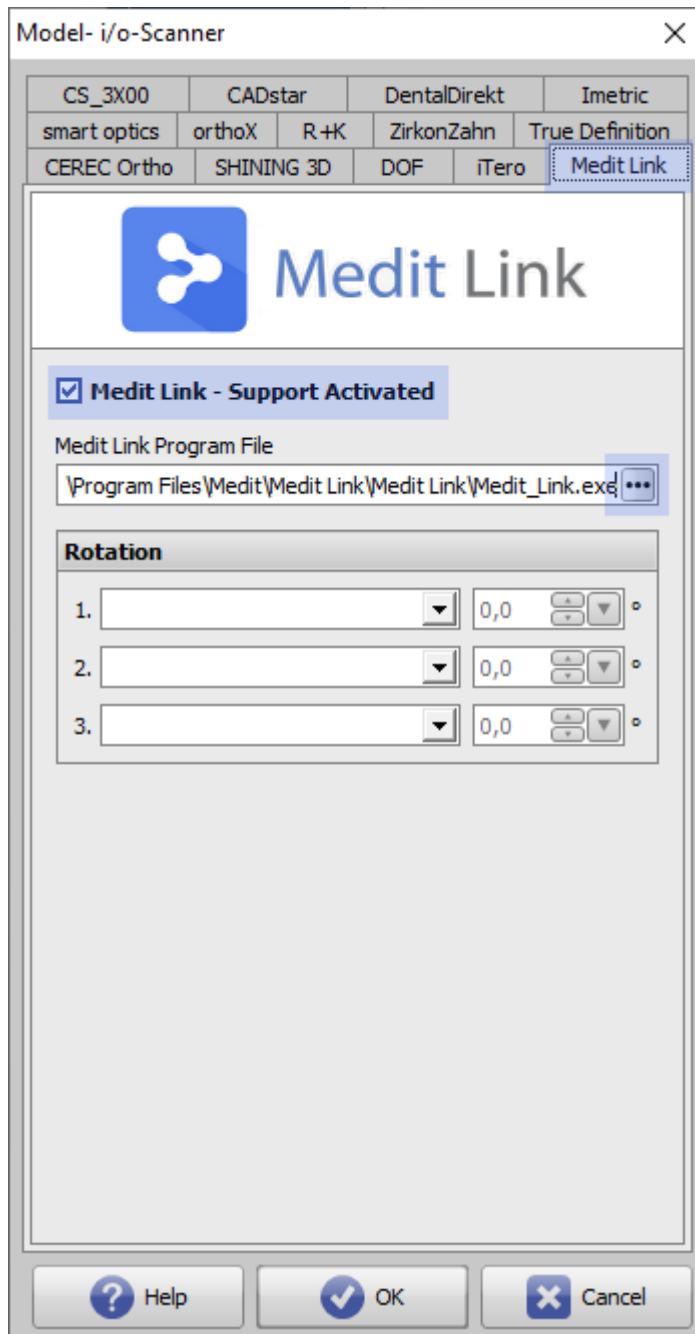


Transfer Scan (POST)

1. Open Medit Link on local OnyxCeph^{3™} client computer
2. Select patient | case | scan ①
3. Click icon button [Export to OnyxCeph]


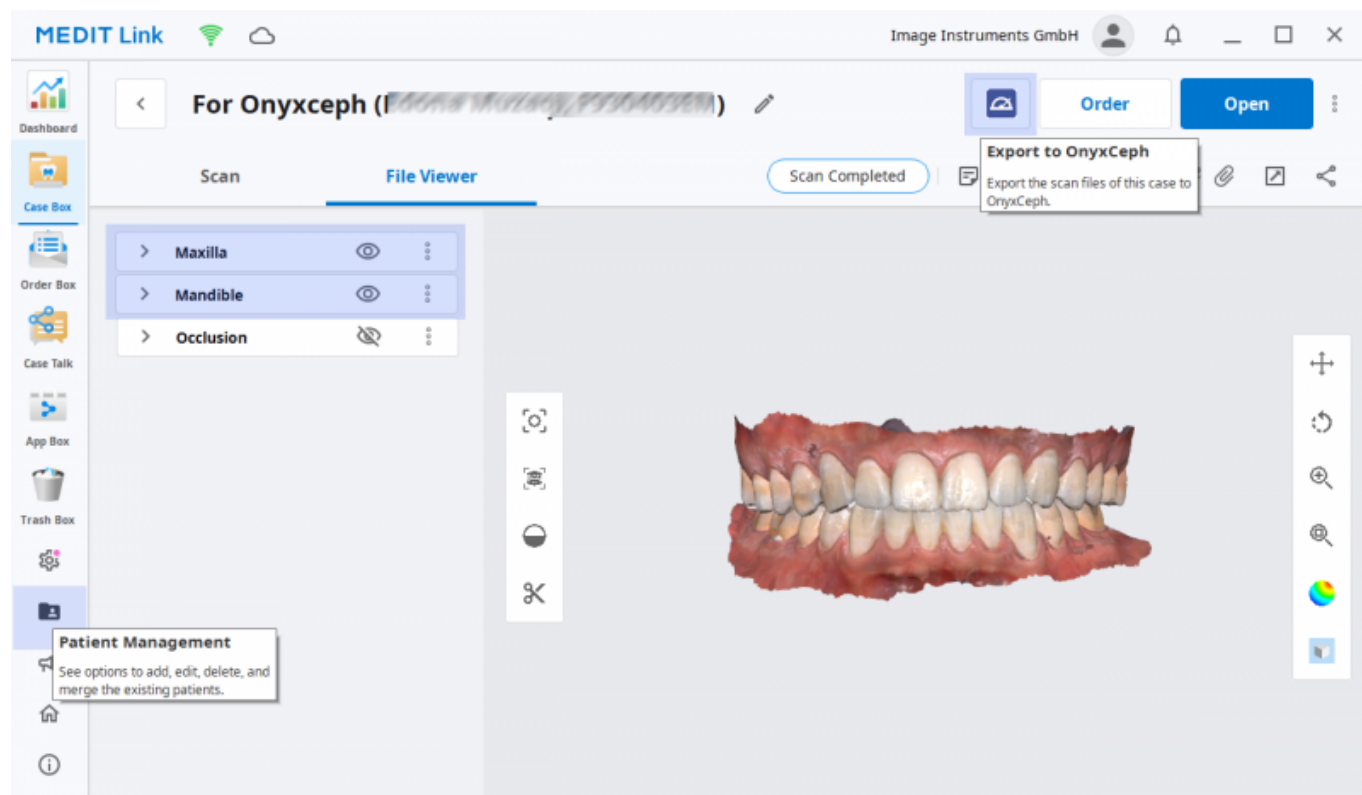


① Patient must be created in the Medit Link patient administration with all master data! It is recommended to transfer the master data from OnyxCeph^{3™} or, if supported, from the PMS to ensure a unique assignment of the scans to the patient.

4. If the scan is not associated with the current login, a warning will appear:

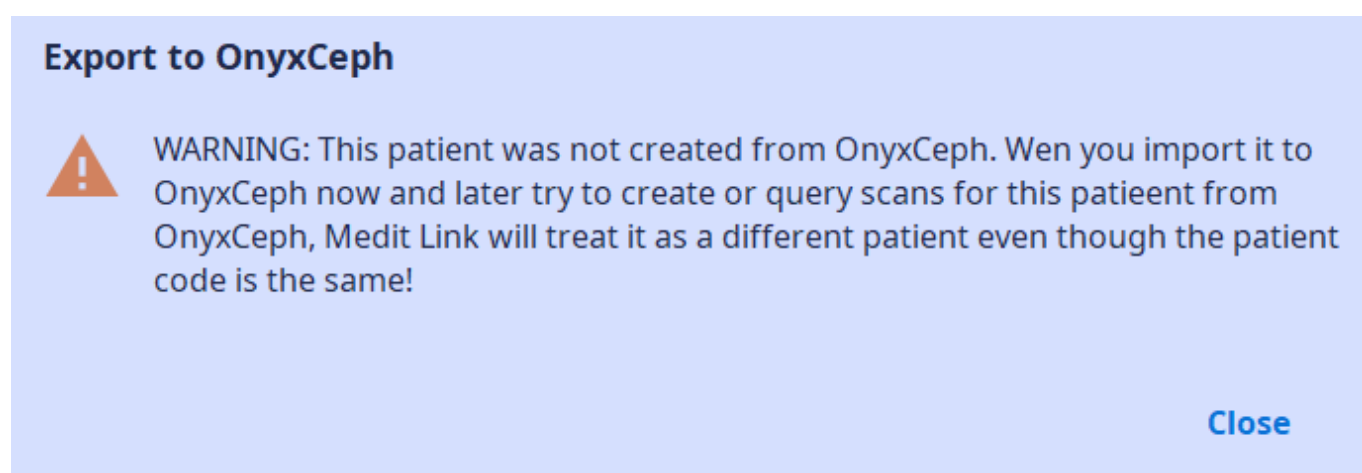


5. If this warning is closed, the scan data is transferred to module [Add image 3D](#) and can be classified there, aligned to the patient and further processed according to the known workflow.

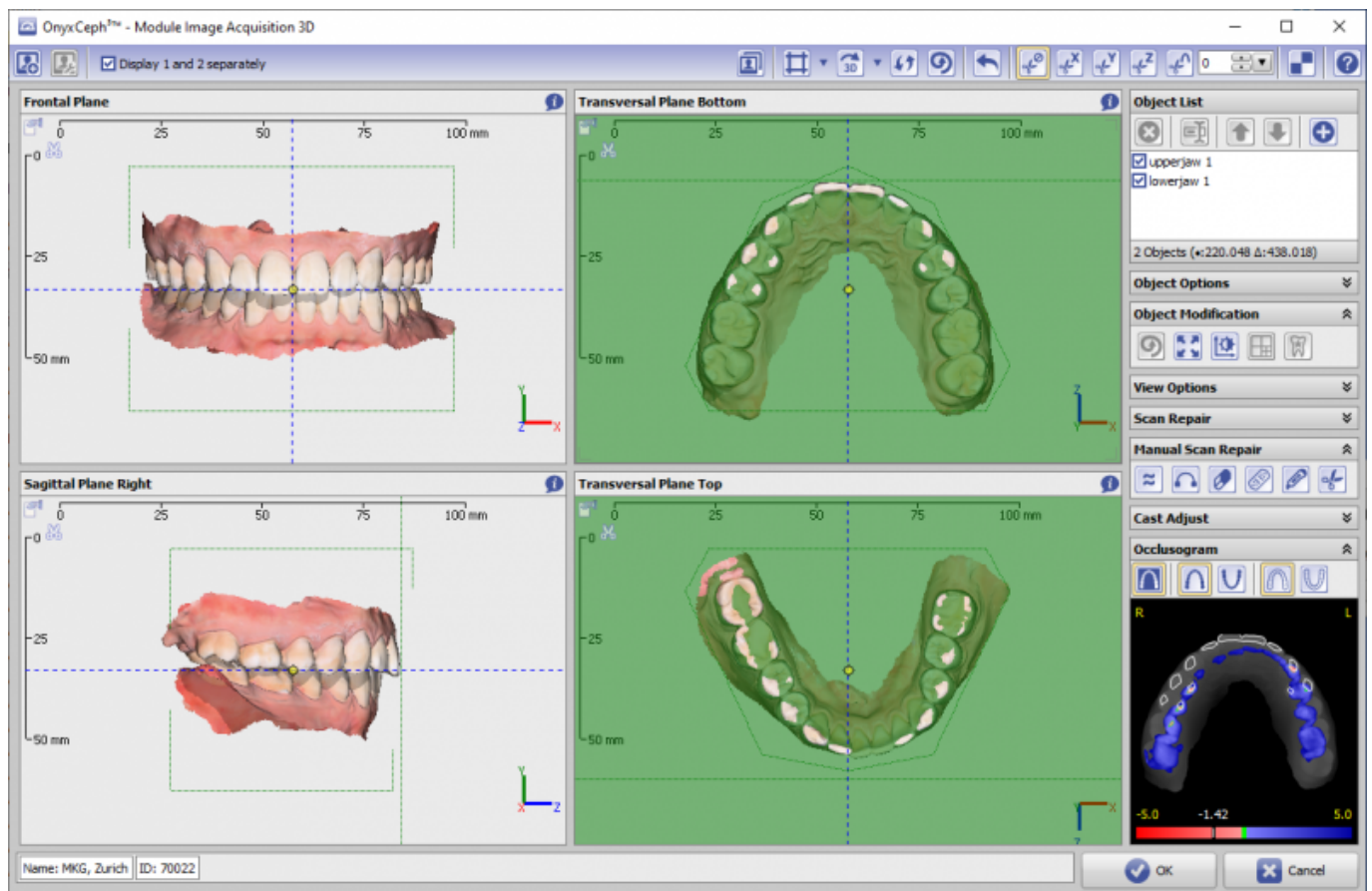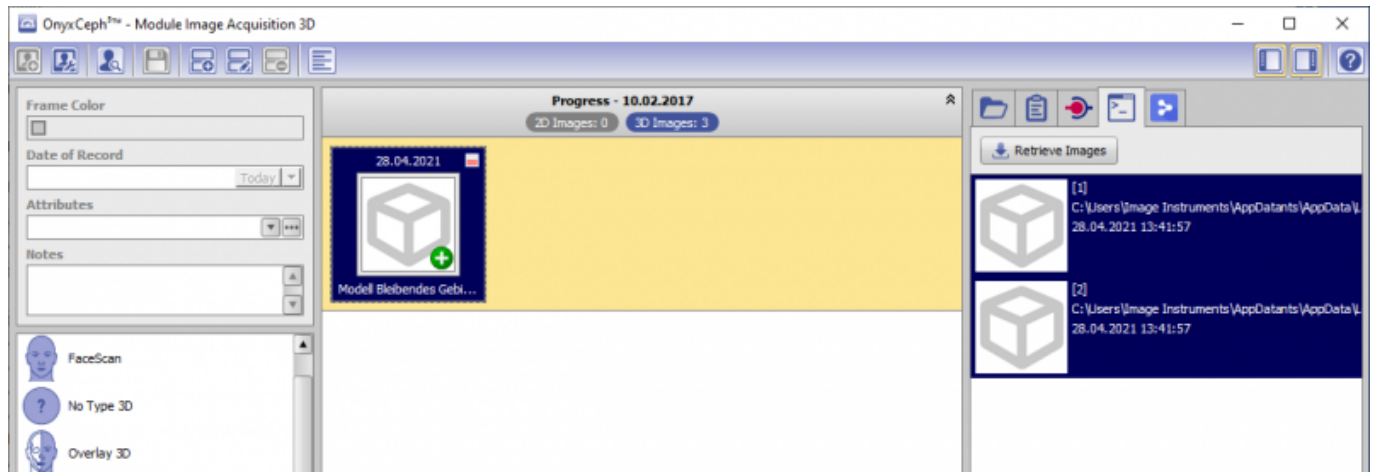


Request Scan (GET)


1. Select patient record in OnyxCeph³™
2. Open module [Bild hinzufügen 3D](#)
3. In the source panel, select tab [Medit Link]
4. Click button [Image Acquisition...] klicken

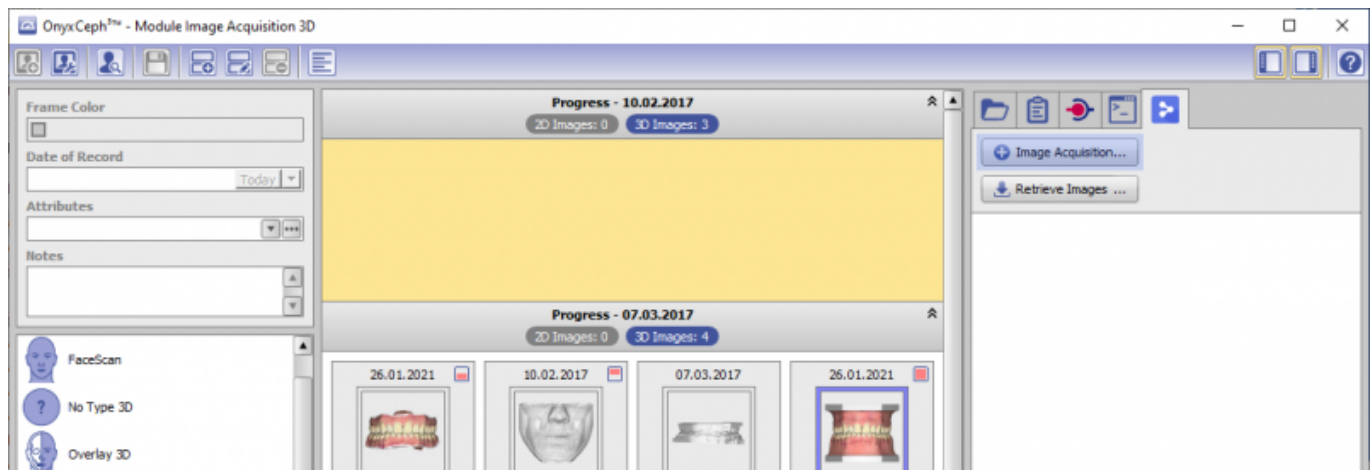


5. In Medit Link, acquire or select scan data

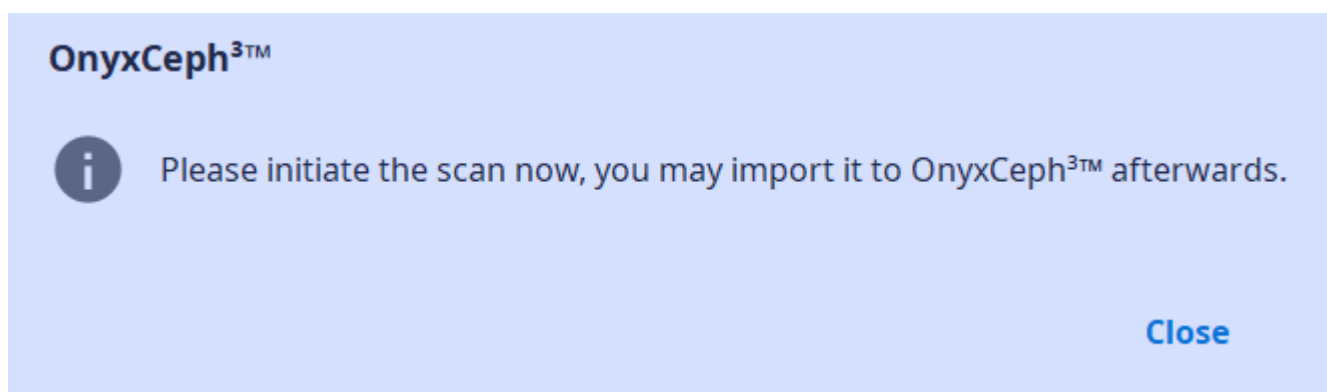

6. List by button [Retrieve Images...] and import, classify, align and further process the scan data according to the known workflow.

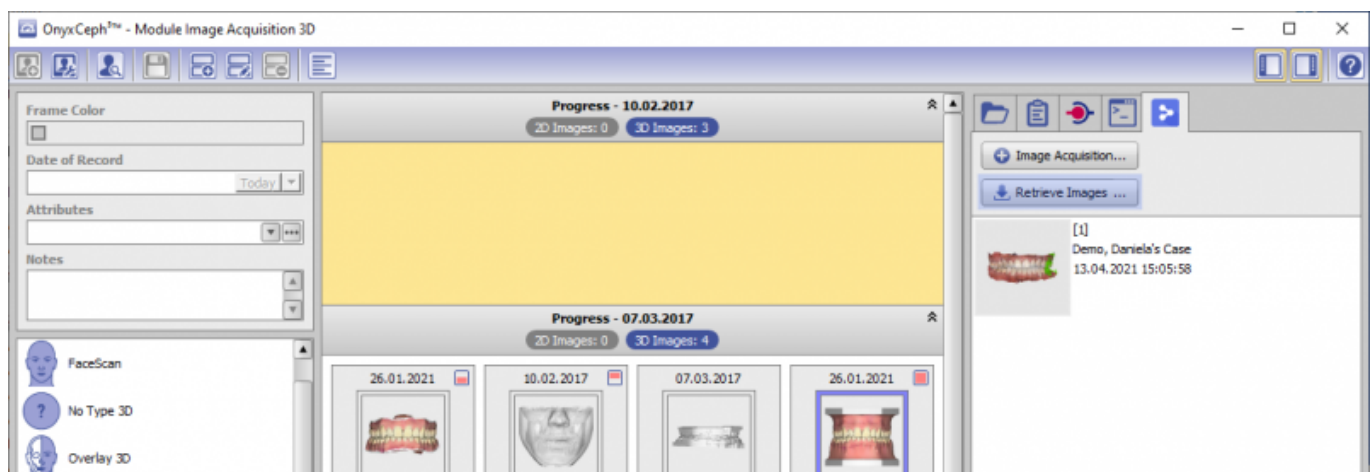

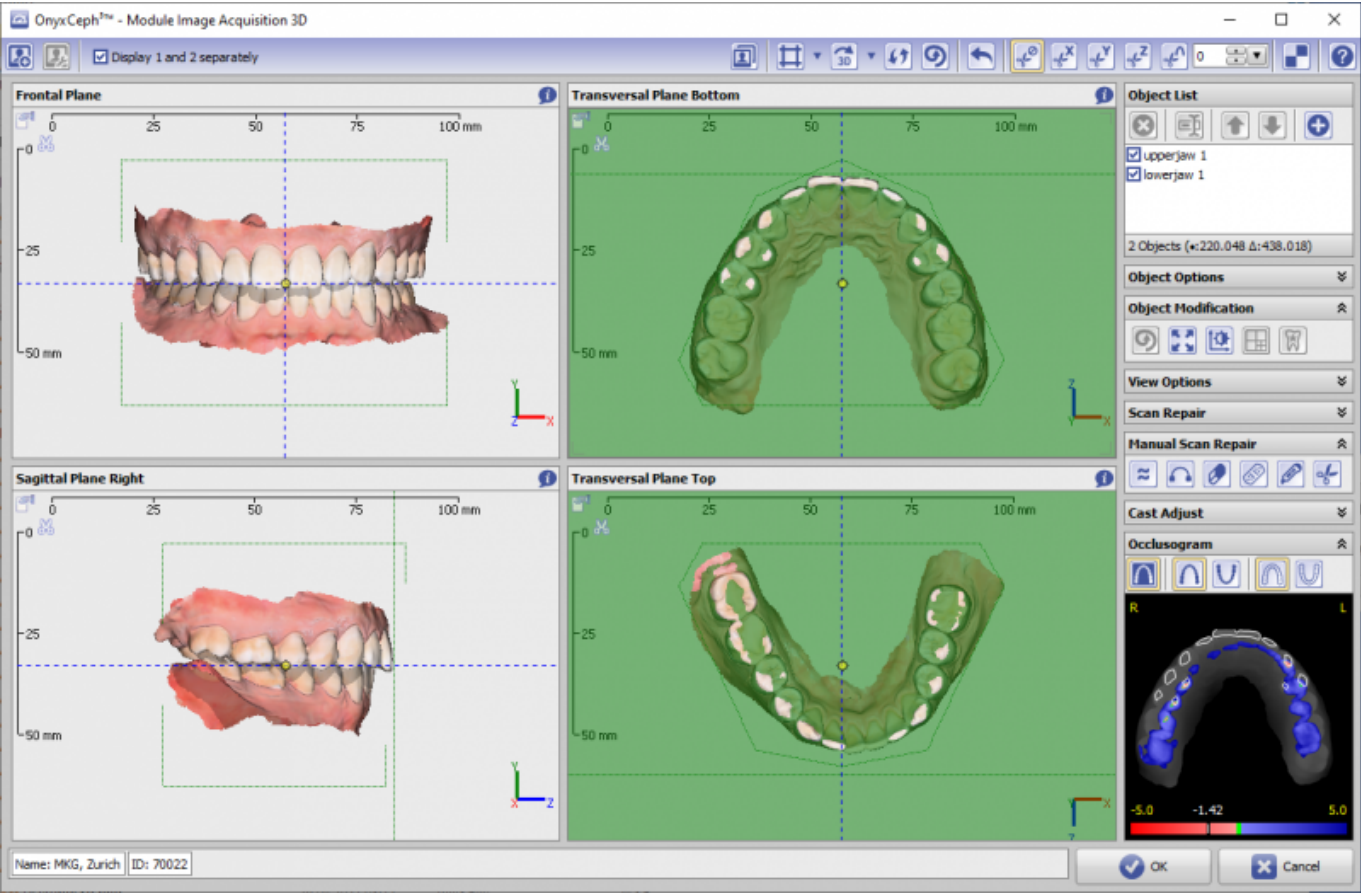

Links



From:

<http://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

<http://onyxwiki.net/doku.php?id=en:meditimport&rev=1655442705>

Last update: **2022/06/17 07:11**

