


Module Combine 2D

Module Combine 2D can combine two images with different but overlapping content. The second image behind the first image is shown completely or in a self-defined region of interest using common reference points with selectable transparency. In addition, the tracing of one of the two images can be overlaid.


First Steps

1. Open module with 2 images
2. Select or digitize correspondin landmarks
3. Define region of interest
4. Switch to tab |Combine|
5. Optionally, modify superimposition
6. Save result as project and/or as finding

Module user interface

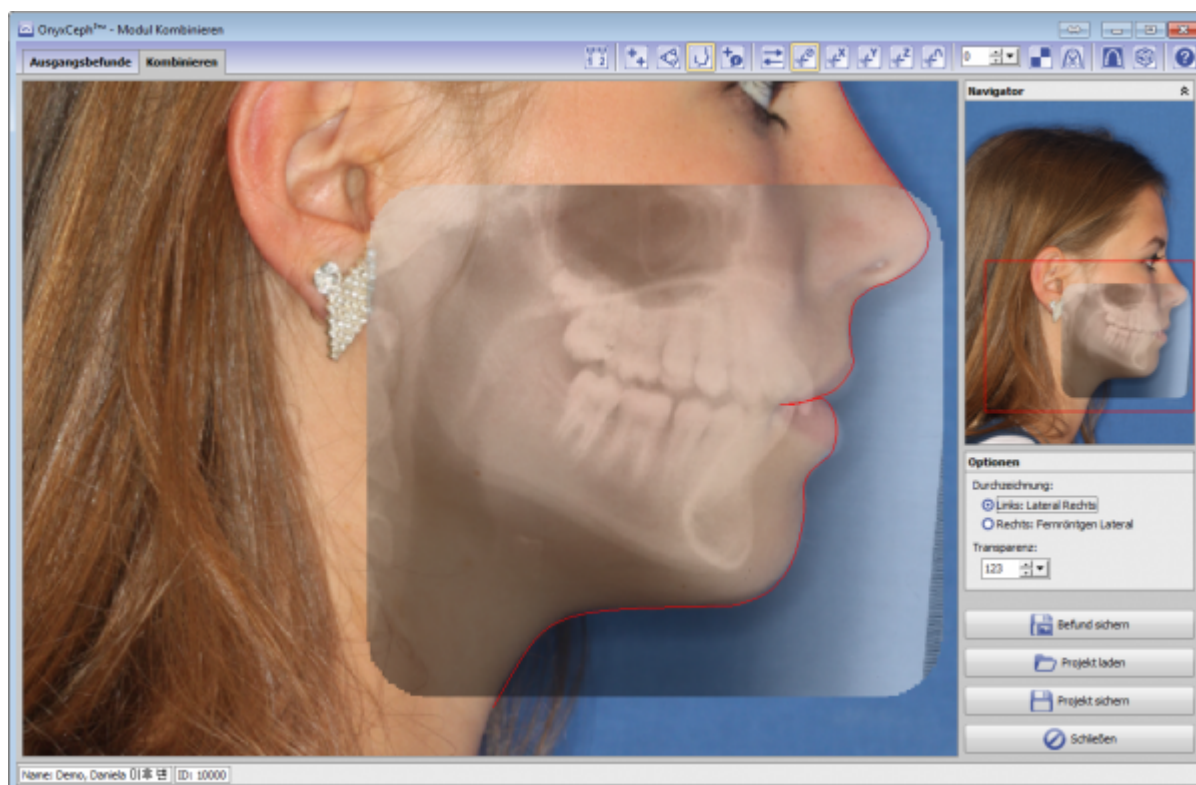

Links

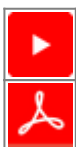

From:

<https://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:

https://onyxwiki.net/doku.php?id=en:combine_2d&rev=1584733993

Last update: **2020/03/20 20:53**

