


System Options - Tab |Printouts|

On system options tab |Printouts|, several settings regarding print templates and their processing and printing can be done.





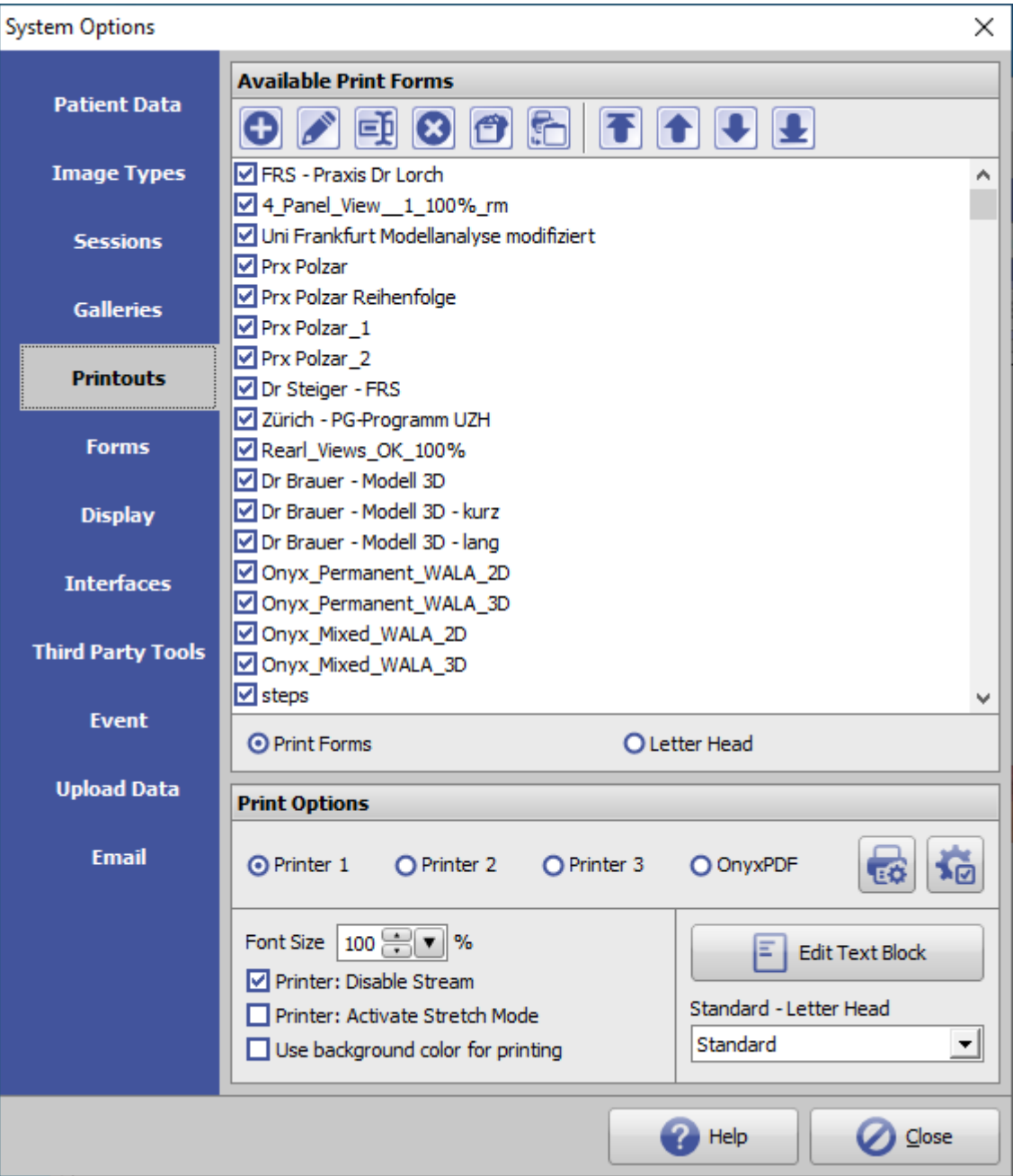




[Click on tab to switch to related topic]

Panel Available Print Forms

► List view of all available templates of type

- Print Forms
- Letter Head

► *Buttons* for lis entry modification

 Up	Moves selected list entry one position up - if [CTRL] is pressed: to first position
 Down	Moves selected list entry one position down - if [CTRL] is pressed: to last position
 Add	Opens editor with empty template
 Edit	Opens editor with template of for selected list entry
 Delete	Deletes for selected list entry
 Export	Exports for selected list entry as patch file
 Convert	Converts print template to presentation template

Panel Print Options

- *Radio buttons Printer*

Printer-related presettings can be made here.

- Preset printer type and printer settings for max 3 Windows printers and Onyx internal PDF printer
- Definition of font and line parameters for max 3 Windows printers and Onyx internal PDF printer



Print Options

☒ Printer 1

☐ Printer 2

☐ Printer 3

☐ OnyxPDF



Print

Printer

Name:

Microsoft Print to PDF

Properties...

Status:

Ready

Type:

Microsoft Print To PDF

Where:

PORTPROMPT:

Comment:

☐ Print to file

Print range

☒ All

☐ Pages from: to:

☐ Selection

Copies

Number of copies:

1

1

2

3

☐ Collate

OK

Cancel

View

Input Window



Font

Font

☐ Shadow

☒ Frame

Lines

Style 1

Size

1

Color

☐ Dotted

Help

OK

Cancel

- *Feature list*


In addition, further print options can be preset:

Font Size %

☒ Printer: Disable Stream

☐ Printer: Activate Stretch Mode

☐ Use background color for printing

 Edit Text Block

Standard - Letter Head

Font Size	Multiplier of the preset system font size
Printer: Disable Stream	Reloads the saved printer setting from the INI
Printer: Activate Strech Mode	Not used anymore
User background color for printing	Layout-background colour [not] used for printing

- **Button Edit Text Block**

The text block defined here can be inserted into text fields under category [Patient] in the [print/presentation editor](#).

Editor

Patient Data

```
{pat:first name} {pat:name}
{pat:short_age}
{pat:gender}
{pat:number}
```

Line: 5 Column: 9 Ins

- Pulldown Standard Letter Head

The letterhead selected here will be used for print form "Standard".

From:
<https://onyxwiki.net/> - [OnyxCeph³™ Wiki]

Permanent link:
https://onyxwiki.net/doku.php?id=en:tab_printouts

Last update: **2021/03/25 14:40**

